

HORIZON

12 Flex, 24 Flex









Operator's Manual



TABLE OF CONTENTS

Symbols	1
Caution Warning Statements	2
Model Description	4
Features.....	4
Intended Use	4
Warranty	4
Initial Setup.....	5
Quick Start	6
Settings	7
Quick Adjust Time and Speed.....	7
Adjusting the Brake Setting	7
Changing the Beeper (Audible Alert).....	8
Create New Cycle	8
Display Cycle Count	8
Modifying a Preset	9
Deleting a Cycle	10
Cycle Lock.....	10
Preset Cycles.....	11
Operation	11
Balancing Loads	12
12-Tube Centrifuges	12
24-Tube Centrifuges	12
Care and Preventative Maintenance	13
Cleaning and Disinfection	14
Troubleshooting	15
Calculating the G-Force	16
General Specifications	17
Replacement Parts.....	18
Previous Manual Revisions	18

Symbols

Symbol	Definition	Use
	Caution	Caution to safety hazard. Potential risk of personal injury or damage to the instrument if improperly handled. Consult the manual before proceeding.
	Manufacturer	Manufacturer of record.
	Electrical and electronic products recycling symbol	Recycle only as electronic waste. Do not dispose in normal waste.
	RoHS Compliant	Compliance with RoHS environmental standards.
	CE Mark	Denotes conformity to specific European directives and regulations.
	UK Mark	Denotes conformity to specific UK directives and regulations.
	MET Listing	Denotes conformity to specific safety standards and regulations.
FDA LISTED	FDA Listed	Denotes that the product has been properly listed with the FDA.
	ISO Certification	Denotes conformity to quality standards and quality management systems.

CAUTION WARNING STATEMENTS



This device is intended to be operated by properly trained personnel who have carefully read the operating manual and are familiar with the function of the device. Users should also comply with the specimen receptacle manufacturer's specific instructions for use, in addition to any other protocols established by the testing organization.



WARNING: For the safety of both the operator and service personnel, care should be taken when using this centrifuge if handling substances that are known to be toxic, radioactive or contaminated with pathogenic microorganisms. Use appropriate personal protection equipment (PPE). When Risk Group II materials are used, (as identified in the World Health Organization "Laboratory Bio-Safety Manual"), a Bio-Seal should be employed. In the event that materials of a higher risk group are being used, more than one level of protection must be provided. The use of flammable or explosive materials as well as those materials which have a vigorous chemical reaction is prohibited.



Unplug the centrifuge before cleaning or performing maintenance.



WARNING: Inspect centrifuge for cracks or physical damage to cabinet, lid, rotor, or tube holders. Damage may result in unsafe operation. Discontinue use until repairs have been performed.



The use of flammable or explosive materials as well as those materials which have a vigorous chemical reaction is prohibited.



For your safety and durability of the machine, never transport or store centrifuge with tube holders inside the machine.



WARNING: "Universal precautions"¹ should be followed in handling all items contaminated with blood or other bodily fluids.



This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with this operator manual, may cause interference to radio communications.



Operation of this equipment in a residential area may cause interference, in which case the user will be required to correct the interference at his own expense.



Operation of this equipment in a manner not specified by the manufacturer may impair the protection provided by the equipment.



Electrical Safety protection is provided by properly connecting the centrifuge to earth ground. Use only the manufacturer provided line cord and ensure that it is connected to a properly grounded power receptacle. Failure to do so will result in an electrical hazard.



WARNING: Do not make modifications to or remove any hardware from rotor without prior authorization from Drucker Diagnostics.



WARNING: Only use Drucker Diagnostics components in this centrifuge.



Due to the lack of the possibility of human exposure, all Drucker centrifuges and accessories sold by Drucker Diagnostics, Inc. are compliant without any special labeling required by the California Safe Drinking Water and Toxic Enforcement Act (Proposition 65).

1 Recommendations for Prevention of HIV Transmission in Health Care Settings. MMWR 1987; 36 (Supplement #2S)

MODEL DESCRIPTION

HORIZON is a versatile line of centrifuges designed with 3 settings to process Chemistry, Coag or Platelet Poor Plasma, and Urine specimens in the same unit. The maximum g-force of 2,000 xg makes HORIZON compatible with most brands of tubes. Cycle settings can be changed to accommodate custom settings.

This general-purpose laboratory centrifuge may also be used to spin approved containers with biologics, chemicals (non-flammable, non-explosive, non-volatile, and non-highly reactive), and environmental samples.

FEATURES

- The first three (3) cycles are conveniently pre-set and labeled for your lab's most common applications. Use the default cycles or settings can be customized.
- A time and speed/g-force can be quickly entered for a single use cycle. The cycle will not be retained in memory.
- If desired, the control panel can be temporarily locked on one cycle, ideal for standardization to a single spin.
- A Preset Lock can be turned on to prevent changes from being made accidentally to programmed cycles.
- Up to 10 cycles can be programmed for time, speed, and braking and labeled with a custom name. Cycles can be programmed by g-force (RCF) or speed (RPM) to facilitate matching validated cycles and tube manufacturers' IFUs.
- A digital cycle counter tracks the number of cycles the centrifuge has run.
- Lid lighting indicates the centrifuge's status (ready, running, done), informing the operator when tubes are ready for the analyzer and preventing tubes from being left in the centrifuge longer than necessary (patent pending).
- A traditional audible alert indicates the completion of the cycle. The audible alert can be muted.
- Cool-Flow design prevents overheating of samples by using ambient air to keep specimens at room temperature.
- The tube holders are fiber reinforced for high strength, durability, and years of trouble-free use.
- A clear lid permits safe observation of samples and optical calibration of speed.
- The lid safety system prevents the centrifuge from operating unless the lid is closed and latched.
- The lid safety system only allows entry into the centrifuge after the rotor has completely stopped.
- The high-power brushless motor provides years of operation with no routine maintenance.

INTENDED USE

General purpose laboratory centrifuge, intended for the density-based separation of fluids through centripetal acceleration.

WARRANTY

Drucker Diagnostics warrants that this centrifuge is free from defects in workmanship and parts for 2 years.

INITIAL SETUP

- Unpack and verify that all the following are included:
 - Centrifuge
 - Power cord
 - Tube holders
 - Quick Start Guide
- Setup the centrifuge on flat and level surface. A bench top clearance height of 21" (54 cm) is required to open the lid.
- The centrifuge should have 6" (15 cm) of clear space around the centrifuge. Proper ventilation is necessary to prevent the overheating of samples as well as premature failure of the centrifuge. Choose an area which allows unencumbered air flow, and where the temperature remains between 16°C and 32°C.
- No hazardous material shall be permitted in the clearance envelope during operation.
- The operator time within the envelope shall be limited to the time necessary for loading, unloading, and centrifuge operation only.
- Plug the line cord into the centrifuge.
- Plug the line cord into an approved electrical outlet.




BE SURE THE ELECTRICAL OUTLET IS ALWAYS ACCESSIBLE AS THE LINE CORD IS THE MEANS OF EMERGENCY DISCONNECTION!

QUICK START

The top screen display alternates between the name of the currently selected speed. The bottom screen displays the time setting. For convenience, the first few cycles are preset as shown and can be reprogrammed to your validated cycles:

- (1) Chemistry This setting is factory preset for Chemistry tubes
- (2) Coag (PPP) This setting is factory preset for Coagulation or Platelet Poor Plasma (PPP)
- (3) Urine This setting is factory preset for urine tubes

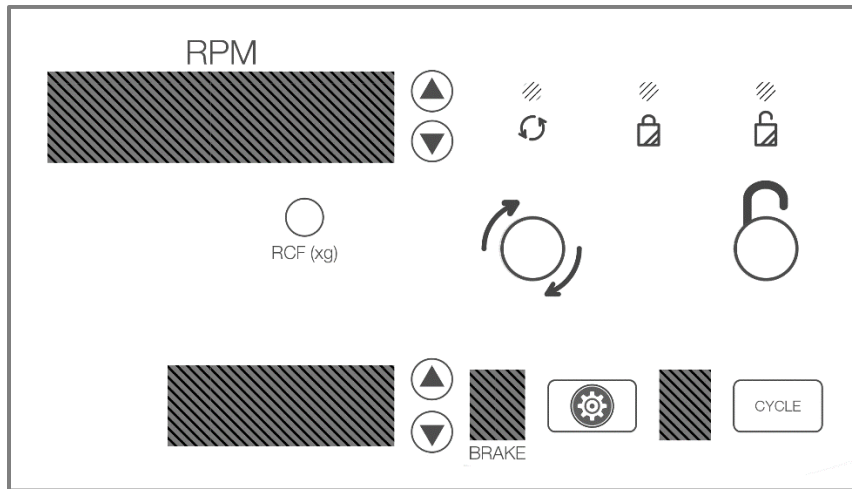




	Start	Begins running the cycle displayed on the screen. The lid must be closed.
	Unlock	Allows access into the rotor chamber by engaging the unlocking mechanism. Entry is only possible when the rotor is stopped.
	Stop	Pressing the UNLOCK button during operation will terminate the run and unlock the lid after the rotor has come to a stop.
	Cycle Selection	Press the CYCLE button to select the desired saved cycle.




SETTINGS

QUICK ADJUST TIME AND SPEED




Change time, speed (RPM) or g-force (RCF) for a single cycle.

	<p>Setting Speed</p>	<p>To change the speed (RPM) shown on the top display, use the up and down buttons next to that screen. The CYCLE number is replaced with a "--" in the display, and the top screen displays the speed.</p>
	<p>Setting by G-Force</p>	<p>Press and hold the RCF (xg) button while changing the displayed setting on the top screen, using the up and down buttons next to it. The RPM will automatically adjust.</p>
	<p>Setting Time</p>	<p>Press the up and down buttons next to the TIME display.</p>



ADJUSTING THE BRAKE SETTING

	<p>Enter the Advanced Menu</p>	<p>Press the GEAR button to enter the advanced menu.</p>
	<p>Change Brake Values</p>	<p>While in the advanced menu, navigate to "Brake". Use the UP and DOWN buttons next to the TIME screen to adjust brake to desired value between 0 (no brake applied) and 9 (maximum braking force applied).</p>
	<p>Exit the Menu</p>	<p>Press the GEAR button.</p>


CHANGING THE BEEPER (AUDIBLE ALERT)

	<p>Enter the Advanced Menu</p>	<p>Press the GEAR button to enter the advanced menu.</p>
	<p>Turn Beeper On or Off</p>	<p>While in the advanced menu, navigate to “Beeper”. Switch ON or OFF with the UP and DOWN buttons next to the TIME display. This setting will apply to all cycles.</p>
	<p>Exit the Menu</p>	<p>Press the GEAR button.</p>

CREATE NEW CYCLE









	<p>Change Settings</p>	<p>Refer to previous table (Quick Adjust Time and Speed) to change speed and time to desired values.</p>
	<p>Save Cycle</p>	<p>Hold the CYCLE button until you hear a double beep.</p>

DISPLAY CYCLE COUNT




	<p>Display Cycle Count</p>	<p>With the lid open and the unit powered, press and hold the START button. The cycle count will be displayed until the START button is released.</p>
---	-----------------------------------	---

MODIFYING A PRESET

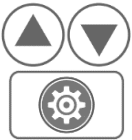Save up to 10 custom cycles. The top screen alternates between cycle name and speed.




	Access the Menu	With desired cycle selected, access the menu.
	Setting by G-Force (Recommended)	Press and hold the RCF (xg) button while changing the setting, using the UP and DOWN buttons next to the display. The RPM will automatically adjust.
	Setting Speed (Alternate)	To change the speed (RPM), use the UP and DOWN buttons next to the display. The g-force will adjust automatically and can be verified by pressing the RCF button.
	Setting Time	Press the UP and DOWN buttons next to the TIME display.
	Change Brake Values	While in the advanced menu, navigate to "Brake". Use the UP and DOWN buttons next to the TIME screen to turn brake on and off.
	Turn Beeper On or Off	While in the advanced menu, navigate to "Beeper". Switch ON or OFF with the UP and DOWN buttons next to the TIME display. This setting will apply to all cycles.
	Naming the Cycle	While in the advanced menu, navigate to the cycle name with the UP and DOWN arrows. Press the START button. The * indicates the space selected. Use the UP and DOWN buttons to change characters, then move to the next space with the right arrow >. Press the GEAR button to return to the main programming menu.
	Save and Exit Settings Mode	Press the GEAR button, followed by the CYCLE button to exit the menu.

DELETING A CYCLE

	<p>Enter the Advanced Menu</p>	<p>With the desired cycle selected, access the menu and enter the advanced menu.</p>
	<p>Navigate to Delete</p>	<p>Using the UP and DOWN buttons, navigate to DELETE. Exit the menu. WARNING: CYCLE WILL BE DELETED IF MENU IS EXITED WITH DELETE SELECTED</p>
	<p>Confirm Deletion</p>	<p>Press the CYCLE button to Delete the cycle</p>

CYCLE LOCK

To ensure repeatability, the centrifuge can be locked either on one cycle (Single Cycle Lock) or restricted to the saved cycles (Presets Only Mode).

	<p>Presets Only Mode</p>	<p>Allows selection of any saved cycle and prevents changing the parameters of saved cycles</p>	<p>With the lid open, press and hold the UNLOCK button to enter Preset Lock. One beep will confirm that cycle selection is now locked. Select desired cycle. NOTE: If Single Cycle Lock is set, it must be canceled before Presets Only Mode can be set.</p>
	<p>Single Cycle Lock</p>	<p>Prevents making changes to the single selected cycle parameter</p>	<p>With the lid open, select the cycle you wish to lock, press, and hold the UNLOCK button. Two beeps will confirm that cycle selection is locked.</p>
	<p>Cancel Lock</p>		<p>Hold the UNLOCK button. Three beeps will confirm that the cycle selection is now unlocked.</p>

PRESET CYCLES

		(1) Chemistry	(2) Coag (PPP)	(3) Urine
HORIZON 12 Flex	RPM	3,500	3,200	1,800
	Time	10	15	5
	G-Force	1,800 xg	1,500 xg	500 xg
	Brake	6	3	3
HORIZON 24 Flex	RPM	3,300	3,000	1,700
	Time	10	15	5
	G-Force	1,800 xg	1,500 xg	500 xg
	Brake	6	3	3

OPERATION

- Place the tubes into the tube holders. Be sure to follow the rules for balanced loads as listed in next section.
- Close the lid and turn the lid knob clockwise to its complete stop position.
- The digital screen shows the currently selected cycle. To select another cycle, press the UP or DOWN button in succession until the desired cycle is selected.
- Pushing the START button on the control panel will start the spin cycle.
- When the cycle is completed, the rotor will slow to a complete stop and the lid light will flash.
- The unlocking mechanism will engage for 60 seconds allowing entry into the rotor chamber. To unlock after more than 60 seconds have elapsed, press the UNLOCK button. The lid will unlock for another 15 seconds.
- Turn the lid knob counterclockwise and open the lid. The lid light will turn off.
- You may now safely remove the samples.


BALANCING LOADS



Your centrifuge must contain a balanced load to work properly. Spinning balanced loads will extend the life of the centrifuge and produce better results. Use the following rules when loading the rotor. If an odd number of samples is to be spun, fill a tube with water to match the weight of the unpaired sample and place it across from this sample.

Opposing tube holders must be equally loaded or empty or loaded with equally weighted samples. When loading only 3 tubes, they must be of equal weight.


12-TUBE CENTRIFUGES

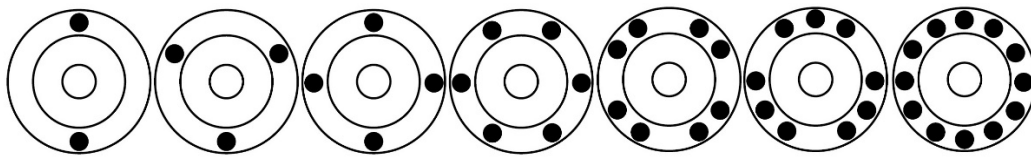


24-TUBE CENTRIFUGES

Buckets can be placed around the rotor in any of the rotor loading configurations shown below. Each bucket must be loaded symmetrically with tubes.


Balanced Rotor Loading

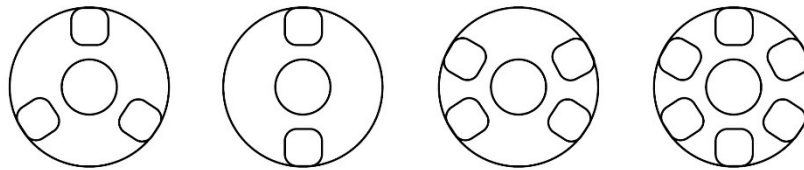

Proper Bucket Loading

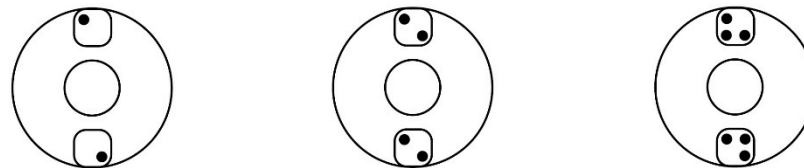

15 mL / 125 mm tubes may be spun in the standard HORIZON 24 Flex buckets, but must be loaded as follows:

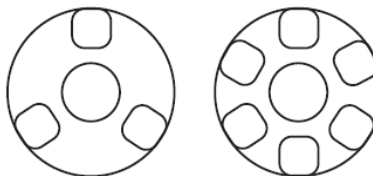

When three buckets are being used, each bucket must be loaded symmetrically in any configuration below:



When all six buckets are being used, the rotor must be loaded in the following configuration **only**. Opposing buckets do not need to be filled but must be loaded in exactly the same way.



CARE AND PREVENTATIVE MAINTENANCE

With proper care and maintenance, your centrifuge will provide years of laboratory service. For proper care, the following steps should be taken:

- **Always Spin Balanced Loads:** Make certain that you are always spinning a balanced load, as shown in the previous section. These centrifuges have a unique counter balanced motor mounting design which produces excellent vibration dampening. However, out-of-balance loads may break glass test tubes and may produce unsatisfactory separation results. Proper load balancing will improve sample separation and extend the life of the centrifuge.
- **Motor and Electrical Maintenance:** The highest quality electrical components have been selected for the centrifuges and should not need maintenance or servicing for the life of the centrifuge.
- **Tube Holder Replacement:** It is recommended that the tube holders be replaced after 24 months of use. Inspect tube holders regularly for cracks. If cracks are discovered, replace immediately.
- **Remove Accessories Before Moving:** All tube holders, samples, and caps must be removed from the rotor chamber before transporting or storing the centrifuge to prevent damage and injury.

CLEANING AND DISINFECTION

To prolong the life of the centrifuge, cleaning and disinfection is recommended every six months or whenever there is a spillage or tube breakage. Contaminants must be removed immediately or corrosion and premature degradation of components can occur. Before using any cleaning or decontamination methods other than those recommended by the manufacturer, users should verify with the manufacturer that the proposed method will not damage the equipment.



Cleaning and Decontamination may be necessary as a safeguard before laboratory centrifuges, rotors, and any accessories are maintained, repaired, or transferred.

- Unplug the centrifuge before cleaning.
- Use appropriate personal protective equipment (PPE).
- Apply Cleaning solutions with a dampened towel or cloth **ONLY**. Do not spray or pour cleaning solution directly onto or into the centrifuge. Do not saturate or submerge the centrifuge in water or other cleaning solutions as this will cause damage, create a safety risk, and void the warranty.
- **ONLY** isopropyl alcohol or a 10% (5500 PPM) bleach solution should be used to disinfect the centrifuge and its accessories.
- All surfaces must be dried immediately after cleaning and disinfecting.




TBQ GERMICIDAL PRODUCTS ARE NOT RECOMMENDED AS THEY MAY CAUSE DAMAGE TO THE CENTRIFUGE. REFRAIN FROM USING TO PREVENT VOIDING THE WARRANTY.

- Fully/partially halogenated hydrocarbons, ketones, esters, ethers, benzyls, ethyl benzenes, and all other chemicals not prescribed by the manufacturer shall not be used as they may cause damage to the rotor chamber, rotor, tube holders, accessories and centrifuge exterior and void the warranty.

TROUBLESHOOTING

NOTE: The latch must be turned completely clockwise to its stop position for the centrifuge to operate.

<p>The centrifuge does not run</p>	<ul style="list-style-type: none"> ○ Verify that the centrifuge is powered. One of the LED lights should be on. ○ If “Lid not closed” message is displayed, make sure the lid latch is turned completely clockwise to its stop position. ○ If the centrifuge still does not run, contact Customer Service.
<p>The rotor does not spin freely</p>	<ul style="list-style-type: none"> ○ Make sure nothing has fallen into the rotor chamber, following the procedure above. ○ If nothing obstructs the rotor, the rotor may be damaged. Contact Customer Service for further assistance.
<p>The centrifuge makes a rattling noise when running</p>	<ul style="list-style-type: none"> ○ Stop the centrifuge. Open the lid. ○ Wearing PPE, remove tubes and tube holders/buckets and look for fallen objects or debris. Carefully reach inside the rotor chamber with a tool to remove them. ○ Inspect the rotor, tube holders or buckets for damage. ○ If the tube holders or buckets have any damage, even slight, safely dispose of them and replace them. ○ If the rotor appears damaged, contact Customer Service for further assistance.
<p>Excessive noise or vibration when the centrifuge is running</p>	<ul style="list-style-type: none"> ○ Verify that all four centrifuge feet are properly seated on a flat surface. ○ Ensure that the load is balanced according to instructions in the “Balancing Loads” section of this manual. ○ Make sure that nothing has fallen into the rotor chamber.
<p>“Abort” is displayed on the top screen</p>	<p>The centrifugation cycle has been interrupted.</p>
<p>The centrifuge stops and beeps continuously</p>	<p>The load is not balanced. Press the UNLOCK button, open the lid, and balance the load as recommended elsewhere in this manual.</p>
<p>The centrifuge is stuck on one of the settings</p>	<p>Cycle selection is locked. Press the UNLOCK button for 5 seconds.</p>
<p>Only a few cycles can be accessed</p>	<p>The Preset Lock is active. To deactivate it, press the UNLOCK button for 5 seconds, until you hear 2 beeps, then again until the next 2 beeps. All cycles can now be accessed and/or amended.</p>
<p>The cycle time and speed are not set to the desired value</p>	<p>Check the setting by following the instructions in the section on Setting or Modifying a Saved Cycle. If the preset is not the desired length, follow the procedure on the same page to change the preset time.</p>

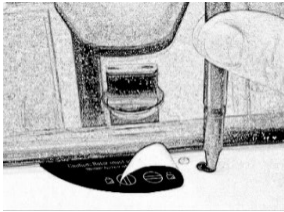Cycle parameters cannot be changed	<ul style="list-style-type: none"> ○ If cycle selection is locked on one cycle, press the UNLOCK button for 5 seconds. Then, press the GEAR button and follow the instructions elsewhere in this manual. ○ If different saved cycles can be selected but not modified, the centrifuge is in Preset Lock mode. Press the UNLOCK button for 5 seconds until two beeps are heard, then again until the next two beeps. You should now be able to change cycle parameters.
The centrifuge does not unlock after a run is completed	<ul style="list-style-type: none"> ○ Wait until the rotor has come to a complete stop. If the lid knob still cannot be rotated, press the UNLOCK button and try again. <ul style="list-style-type: none"> ○ If no LED light is on, the unit is not powered and the lid will not unlock by conventional means. Remove the latch label and use a pen to manually disengage the locking mechanism. Pull the mechanism towards the control panel and then unlatch and open the lid. ○ If the unit is damaged, contact Customer Service for assistance. 
The lid does not open	<ul style="list-style-type: none"> ○ Ensure that the lid knob is turned fully counterclockwise. ○ If the knob cannot be turned counterclockwise, turn it fully clockwise, press UNLOCK, and turn counterclockwise. ○ If the lid remains locked after this and will not unlock, the electronics may have been damaged. Contact customer service for assistance.
Clicking noise during braking gets loud	<ul style="list-style-type: none"> ○ Make sure that the screw in the center of the rotor is tight.
Lid does not stay up	<ul style="list-style-type: none"> ○ Tighten the center screw on the lid hinge.

CALCULATING THE G-FORCE

The I.F.U.s of tube manufacturers recommend cycles at a minimum G-Force, which can be calculated if you know the RPM and the radius. Use the formula below or go to www.druckerdiagnostics.com/g-force-calculator/.

In Centimeters:
 $RCF \text{ or G-force} = 0.00001118 \times \text{Rotor Radius (cm)} \times (\text{RPM})^2$

In Inches:
 $RCF \text{ or G-force} = 0.0000284 \times \text{Rotor Radius (in)} \times (\text{RPM})^2$

	HORIZON 12 Flex	HORIZON 24 Flex
Radius	5.25 in (13.3 cm)	6 in (15.3 cm)

GENERAL SPECIFICATIONS

The rotor and accessories are rated for the maximum rotation frequency shown in the table below.

	HORIZON 12 Flex	HORIZON 24 Flex
Tube Capacity	12 tubes – 3 to 10 mL	24 tubes – 3 to 10 mL, 12 tubes – 15 mL 6 tubes - 50 mL
Dimensions (H x W x D)	13 in x 15 in x 9 in (33 cm x 38 cm x 23 cm)	15 in x 17 in x 9 in (38 cm x 43 cm x 23 cm)
Weight	30 lbs. (13.6 kg)	31 lbs. (14 kg)
Sound Level	56 dB A	59 dB A
Supply Voltage	100 – 240 (+/- 10%)	100 – 240 (+/- 10%)
Supply Frequency	50 - 60 Hz	50 - 60 Hz
Current Consumption	3.6A at 115 VAC 1.8A at 230VAC	3.6A at 115 VAC 1.8A at 230VAC
Centrifuge Motor	½ H.P. Brushless	½ H.P. Brushless
Max g-Force	2,000 xg	2,000 xg
Max Speed	3,700 RPM (+/- 100)	3,400 RPM (+/- 100)
Cycle Time	.5 to 99 minutes (+/- 2%)	.5 to 99 minutes (+/- 2%)
Environmental Conditions		
Set-up Site	Indoor Use Only	Indoor Use Only
Altitude	Up to 2,000m from Sea Level	Up to 2,000m from Sea Level
Ambient Temperature	5 °C to 40 °C	5 °C to 40 °C
Humidity	Maximum relative humidity 80% for temperatures up to 31 °C, decreasing linearly to 50% relative humidity at 40 °C	Maximum relative humidity 80% for temperatures up to 31 °C, decreasing linearly to 50% relative humidity at 40 °C
Overvoltage Category	II	II
Pollution Degree	2	2

REPLACEMENT PARTS

Horizon 12 Flex

Part No.	Description
7724177K	Foot, rubber (Pack of 4)
02-002-1-0044	Lid Tray Assembly
03-1-0006-0008	Motor Assembly
02-006-0-0021	PC Board
7760006	Power cord
03-1-0005-0193	Power Supply, 350W
02-002-1-0041	Lid Assembly
7724071K	Kit, Hinge (Pack of 2)
02-002-1-0057	Seal, lid gasket
03-0-0003-0332	Open/Close Label
7713079K	Kit, 75/100mm Tube Holder, Blue (Pack of 6)
02-006-1-0065K	Kit, Blue LED PCBA Assembly, Blue
03-0-0003-0380	Horizon12 Flex Front Panel Label
00-100-100-010	12/24 Series Soft Button & Spacer Replacement Kit
00-100-100-005	Replacement Grommet & Bushing Kit

Horizon 24 Flex

Part No.	Description
7728052K	Foot, rubber (Pack of 4)
02-006-1-0044	Lid Tray Assembly
03-1-0006-0008	Motor Assembly
02-006-0-0021	PC Board
7760006	Power cord
03-1-0005-0193	Power Supply, 350W
02-002-1-0037	Lid Assembly
7724071K	Kit, Hinge (Pack of 2)
02-002-1-0058	Seal, lid gasket
03-0-0003-0332	Open/Close Label
7713023K	Kit, 4PI Carrier Assembly (Pack of 6)
02-002-1-0065K	Kit, Blue LED PCBA Assembly, Blue
03-0-0003-0379	Horizon 24 flex Front Panel Label
00-100-100-010	12/24 Series Soft Button & Spacer Replacement Kit
00-100-100-005	Replacement Grommet & Bushing Kit

PREVIOUS MANUAL REVISIONS

To access previous manuals, please click the relevant link below:

[Revision F](#)

[Revision G](#)

[Revision H](#)

This operator's manual is part number 03-0-0002-0139 Rev. I

Product Family: HORIZON Series (HORIZON 12 Flex, 24 Flex)

Complies with UL61010-1/CSA C22.2 No. 61010-1 and IEC61010-2-020

Protected by U.S. Patents #6,811,531, # 7,422,554, #D718,463, & #D734,489. Other Patents Pending

FDA LISTED



E112532



INSTRUCTIONS FOR DISPOSAL OF WEEE BY USERS IN THE EUROPEAN UNION



This product must not be disposed of with other waste. Instead, it is the user's responsibility to dispose of their waste equipment by handing it over to a designated collection point for the recycling of waste electrical and electronic equipment. The separate collection and recycling of your waste equipment at the time of disposal will help to conserve natural resources and ensure that it is recycled in a manner that protects human health and the environment. For more information about where you can drop off your waste equipment for recycling, please contact your local city office, waste disposal service, or where you purchased the product.

Designed, built, and supported in the USA



200 SHADY LANE, SUITE 170 – PHILIPSBURG, PA 16866, USA
+1-866-265-1486 (U.S. ONLY) - +1-814-692-7661
CUSTOMERSERVICE@DRUCKERDIAGNOSTICS.COM
DRUCKERDIAGNOSTICS.COM

