

DASH

Flex 24








Operator's Manual



TABLE OF CONTENTS

Symbols.....	1
Model Description	2
Features.....	2
Intended Use	2
Warranty.....	2
Caution and Warning Statements	3
Initial Setup.....	4
Operation.....	4
Quick Start	5
Settings	6
Quick Adjust Time and Speed	6
Adjusting the Brake Setting	6
Changing the Audible Beeper	7
Create New Cycle.....	7
Display Cycle Count	7
Modifying a Preset.....	8
Deleting a Cycle	9
Cycle Lock	9
Preset Cycles.....	9
Balancing Loads	10
Care and Preventative Maintenance	11
Cleaning and Disinfection	11
Troubleshooting	12
General Specifications	14
Calculating the G-Force	14
Replacement Parts.....	15

Symbols

Symbol	Definition	Use
	Caution	Caution to safety hazard. Potential risk of personal injury or damage to the instrument if improperly handled. Consult the manual before proceeding.
	Manufacturer	Manufacturer of record.
	Electrical and electronic products recycling symbol	Recycle only as electronic waste. Do not dispose in normal waste.
	RoHS Compliant	Compliance with RoHS environmental standards.
	CE Mark	Denotes conformity to specific European directives and regulations.
	UK Mark	Denotes conformity to specific UK directives and regulations.
	MET Listing	Denotes conformity to specific safety standards and regulations.
FDA LISTED	FDA Listed	Denotes that the product has been properly listed with the FDA.
ISO13485 certified	ISO Certification	Denotes conformity to quality standards and quality management systems.

MODEL DESCRIPTION

DASH is a versatile line of centrifuges designed for STAT sample processing. Minimize lab error, reduce turnaround time, and streamline training with the simple controls, powerful performance, and unique LED lid lighting system.

As a general purpose laboratory centrifuge, it is designed to also run other approved containers filled with chemicals (non-flammable, non-explosive, non-volatile, and non-highly reactive only), environmental samples, and other non-human body samples.

FEATURES

- The first three (3) cycles are conveniently pre-set and labeled for you lab's most common applications. Use the default cycles or settings can be customized.
- Lid lighting indicates the centrifuge's status (ready, running, done), informing the operator when tubes are ready for the analyzer and preventing tubes from being left in the centrifuge longer than necessary (patent pending).
- Up to 10 cycles can be programmed for time, speed, and braking and labeled with a custom name. Cycles can be programmed by g-force (RCF) rather than speed (RPM) to facilitate matching validated cycles and manufacturers IFUs.
- A time and speed/g-force can be quickly entered for a single use cycle. The cycle will not be retained in memory.
- If desired, the control panel can be temporarily locked on one cycle, ideal for standardization to a single spin.
- A Preset Lock can be turned on to prevent changes from being made accidentally to programmed cycles.
- A traditional audible alert indicates the completion of the cycle. The audible alert can be muted.
- Cool-Flow design prevents overheating of samples by using ambient air to keep specimens at room temperature.
- The tube holders are fiber reinforced for high strength, durability, and years of trouble-free use.
- A clear lid permits safe observation of samples and optical calibration of speed.
- The lid safety system prevents the centrifuge from operating unless the lid is closed and latched.
- The lid safety system only allows entry into the centrifuge after the rotor has completely stopped.
- The high-power brushless motor provides years of operation with no routine maintenance.

INTENDED USE

General purpose laboratory centrifuge, intended for the density-based separation of fluids through centripetal acceleration.

WARRANTY

Drucker Diagnostics warranties that this centrifuge is free from defects in workmanship and parts for 2 years.

CAUTION AND WARNING STATEMENTS

- ⚠ This device is intended to be operated by properly trained personnel who have carefully read the operating manual and are familiar with the function of the device. [Refer to the clinical laboratory method specified by the specimen receptacle manufacturer or established by the medical technology for the products applications.]
- ⚠ For the safety of both the operator and service personnel, care should be taken when using this centrifuge if handling substances that are known to be toxic, radioactive or contaminated with pathogenic microorganisms. Use appropriate personal protection equipment (PPE). When Risk Group II materials are used, (as identified in the World Health Organization “Laboratory Bio-Safety Manual”), a Bio- Seal should be employed. In the event that materials of a higher risk group are being used, more than one level of protection must be provided. The use of flammable or explosive materials as well as those materials which have a vigorous chemical reaction is prohibited.
- ⚠ Unplug the centrifuge before cleaning or performing maintenance.
- ⚠ Inspect centrifuge for cracks or physical damage to cabinet, lid, rotor, or tube holders. Damage may result in unsafe operation. Discontinue use until repairs have been performed.
- ⚠ This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with this operator manual, may cause interference to radio communications.
- ⚠ Operation of this equipment in a residential area may cause interference, in which case the user will be required to correct the interference at his own expense.
- ⚠ Due to the lack of the possibility of human exposure, all Drucker centrifuges and accessories sold by Drucker Diagnostics, Inc. are compliant without any special labeling required by the California Safe Drinking Water and Toxic Enforcement Act (Proposition 65).

INITIAL SETUP

- Unpack and verify that all the following are included:
 - Centrifuge
 - Power cord
 - Tube holders
 - Quick Start Guide
- Setup the centrifuge on flat and level surface. A bench top clearance height of 21" (54 cm) is required to open the lid.
- The centrifuge should have 6" (15 cm) of clear space around the centrifuge. Proper ventilation is necessary to prevent the overheating of samples as well as premature failure of the centrifuge. Choose an area which allows unencumbered air flow, and where the temperature remains between 16°C and 32°C.
- No hazardous material shall be permitted in the clearance envelope during operation.
- The operator time within the envelope shall be limited to the time necessary for loading, unloading, and centrifuge operation only.
- Plug the line cord into the centrifuge.
- Plug the line cord into an approved electrical outlet.
- Turn on the power switch on the back of the centrifuge.

 **BE SURE THE ELECTRICAL OUTLET IS ALWAYS ACCESSIBLE AS THE LINE CORD IS THE MEANS OF EMERGENCY DISCONNECTION!**


OPERATION

- Place the tubes into the tube holders. Be sure to follow the rules for balanced loads as listed on page 10.
- Close the lid and turn the lid knob clockwise to its complete stop position.
- The digital screen shows the currently selected cycle. To select another cycle, press the UP or DOWN button in succession until the desired cycle is selected.
- Pushing the START button on the control panel will start the spin cycle.
- When the cycle is completed, the rotor will slow to a complete stop and the lid light will flash.
- The unlocking mechanism will engage for 60 seconds allowing entry into the rotor chamber. To unlock after more than 60 seconds have elapsed, press the UNLOCK button. The lid will unlock for another 15 seconds.
- Turn the lid knob counterclockwise and open the lid. The lid light will turn off.

You may now safely remove the samples.

QUICK START

The top screen display alternates between the name of the currently selected speed. The bottom screen displays the time setting.

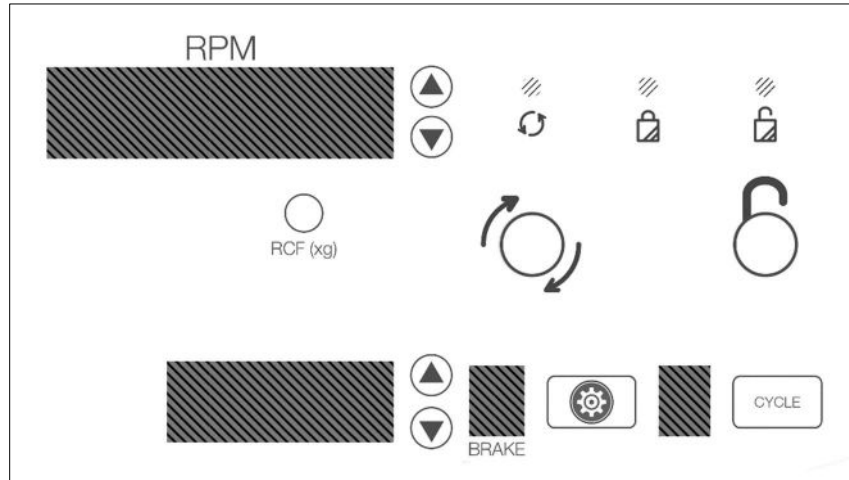




	Start	Begins running the cycle displayed on the screen. The lid must be closed.
	Unlock	Allows access into the rotor chamber by engaging the unlocking mechanism. Entry is only possible when the rotor is stopped.
	Stop	Pressing the UNLOCK button during operation will terminate the run and unlock the lid after the rotor has come to a stop.
	Cycle Selection	Press the CYCLE button to select the desired saved cycle.




SETTINGS

QUICK ADJUST TIME AND SPEED




Change time, speed (RPM) or g-force (RCF) for a single cycle.

	<p>Setting Speed</p>	<p>To change the speed (RPM) shown on the top display, use the up and down buttons next to that screen. The CYCLE number is replaced with a "--" in the display, and the top screen displays the speed.</p>
	<p>Setting by G-Force</p>	<p>Press and hold the RCF (xg) button while changing the displayed setting on the top screen, using the up and down buttons next to it. The RPM will automatically adjust.</p>
	<p>Setting Time</p>	<p>Press the up and down buttons next to the TIME display.</p>



ADJUSTING THE BRAKE SETTING

	<p>Enter the Advanced Menu</p>	<p>Press the GEAR button to enter the advanced menu.</p>
	<p>Change Brake Values</p>	<p>While in the advanced menu, navigate to "Brake". Use the UP and DOWN buttons next to the TIME screen to turn brake on and off.</p>
	<p>Exit the Menu</p>	<p>Press the GEAR button.</p>


CHANGING THE AUDIBLE BEEPER

	<p>Enter the Advanced Menu</p>	<p>Press the GEAR button to enter the advanced menu.</p>
	<p>Turn Beeper On or Off</p>	<p>While in the advanced menu, navigate to “Beeper”. Switch ON or OFF with the UP and DOWN buttons next to the TIME display. This setting will apply to all cycles.</p>
	<p>Exit the Menu</p>	<p>Press the GEAR button.</p>

CREATE NEW CYCLE











	<p>Change Settings</p>	<p>Refer to previous table (Quick Adjust Time and Speed) to change speed and time to desired values.</p>
	<p>Save Cycle</p>	<p>Hold the CYCLE button until you hear a double beep.</p>

DISPLAY CYCLE COUNT




	<p>Display Cycle Count</p>	<p>With the lid open and the unit powered, press and hold the START button. The cycle count will be displayed until the START button is released.</p>
---	-----------------------------------	---

MODIFYING A PRESET

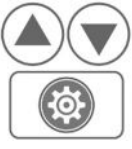Save up to 10 custom cycles. The top screen alternates between cycle name and speed.




	Select Preset Setting	Press the CYCLE button to select the Preset Setting you would like to modify
	Access the Menu	Hold the CYCLE button until you hear a double beep. The cycle number should begin flashing.
	Setting by G-Force (Recommended)	Press and hold the RCF (xg) button while changing the setting, using the UP and DOWN buttons next to the display. The RPM will automatically adjust.
	Setting Speed (Alternate)	To change the speed (RPM), use the UP and DOWN buttons next to the display. The g-force will adjust automatically and can be verified by pressing the RCF button.
	Setting Time	Press the UP and DOWN buttons next to the TIME display.
	Enter the Advanced Menu	Press the GEAR button to enter the advanced menu.
	Change Brake Values	While in the advanced menu, navigate to "Brake". Use the UP and DOWN buttons next to the TIME screen to turn brake on and off.
	Turn Beeper On or Off	While in the advanced menu, navigate to "Beeper". Switch ON or OFF with the UP and DOWN buttons next to the TIME display. This setting will apply to all cycles.
	Naming the Cycle	While in the advanced menu, navigate to the cycle name with the UP and DOWN arrows. Press the START button. The * indicates the space selected. Use the UP and DOWN buttons to change characters, then move to the next space with the right arrow >. Press the GEAR button to return to the main programming menu.
	Save and Exit Settings Mode	Press the GEAR button, followed by the CYCLE button to exit the menu.

DELETING A CYCLE

	Enter the Advanced Menu	With the desired cycle selected, access the menu and enter the advanced menu.
	Navigate to Delete	Using the UP and DOWN buttons, navigate to DELETE. Exit the menu. WARNING: CYCLE WILL BE DELETED IF MENU IS EXITED WITH DELETE SELECTED
	Confirm Deletion	Press the CYCLE button to Delete the cycle

CYCLE LOCK

To ensure repeatability, the centrifuge can be locked either on one cycle (Single Cycle Lock) or restricted to the saved cycles (Preset Lock). The Single Cycle Lock also prevents making changes to the selected cycle parameters. The Preset Lock allows selection of any saved cycle and prevents changing the parameters of saved cycles.

	Enter Preset Lock	Select desired cycle. With lid open, press and hold the UNLOCK button. One beep will confirm that cycle selection is locked.
	Enter Single Cycle	Continue holding the UNLOCK button to enter Single Cycle Lock. Two beeps will confirm that cycle selection is now locked. NOTE: If preset lock is set, it must be canceled before Single Cycle lock can be set
	Cancel Lock	Hold the UNLOCK button. Three beeps will confirm that the cycle selection is now unlocked.

PRESET CYCLES

	Mem-1	Mem-2	Mem-3
RPM	4,300	4,300	3,500
Time	5 min	6 min	7 min
Brake	6	3	3


BALANCING LOADS




Your centrifuge must contain a balanced load to work properly. Spinning balanced loads will extend the life of the centrifuge and produce better results. Use the following rules when loading the rotor. If an odd number of samples is to be spun, fill a tube with water to match the weight of the unpaired sample and place it across from this sample.

*Opposing tube holders must be equally loaded or empty or loaded with equally weighted samples.
When loading only 3 tubes, they must be of equal weight.*

Balanced Rotor Loading

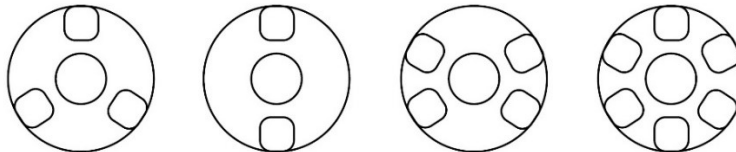


Proper Bucket Loading

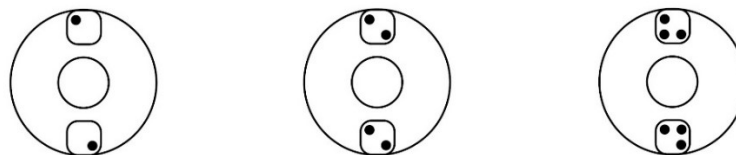

*Buckets can be placed around the rotor in any of the rotor loading configurations shown.
Each bucket must be loaded symmetrically with tubes as above.*

15 mL / 125 mm tubes may be spun in the standard HORIZON 24 Flex buckets, but must be loaded as follows:

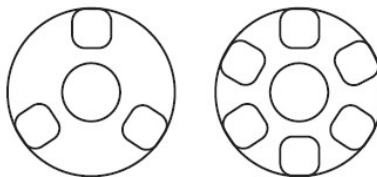

When three buckets are being used, each bucket must be loaded symmetrically in any configuration below:



When all six buckets are being used, the rotor must be loaded in the following configuration **only**. Opposing buckets do not need to be filled but must be loaded in exactly the same way.

CARE AND PREVENTATIVE MAINTENANCE

With proper care and maintenance, your centrifuge will provide years of laboratory service. For proper care, the following steps should be taken:

- **Always Spin Balanced Loads:** Make certain that you are always spinning a balanced load, as shown in the previous section. These centrifuges have a unique counter balanced motor mounting design which produces excellent vibration dampening. However, out-of-balance loads may break glass test tubes and may produce unsatisfactory separation results. Proper load balancing will improve sample separation and extend the life of the centrifuge.
- **Motor and Electrical Maintenance:** The highest quality electrical components have been selected for the centrifuge and should not need maintenance or servicing for the life of the centrifuge.
- **Tube Holder Replacement:** It is recommended that the tube holders be replaced after 24 months of use. Inspect tube holders regularly for cracks. If cracks are discovered, replace immediately.
- **Remove Accessories Before Moving:** All tube holders, samples, and caps must be removed from the rotor chamber before transporting or storing the centrifuge to prevent damage and injury.

CLEANING AND DISINFECTION

To prolong the life of the centrifuge, cleaning and disinfection is recommended every six months or whenever there is a spillage or tube breakage. Contaminants must be removed immediately or corrosion and premature degradation of components can occur. Before using any cleaning or decontamination methods other than those recommended by the manufacturer, users should verify with the manufacturer that the proposed method will not damage the equipment.

- Unplug the centrifuge before cleaning.
- Use appropriate personal protective equipment (PPE).
- Apply cleaning solutions with a towel or cloth. Do not submerge the centrifuge in water or other cleaning solutions as this will cause damage and void the warranty.
- ONLY isopropyl alcohol or a 10% (5500 PPM) bleach solution should be used to disinfect the centrifuge and its accessories.
- All surfaces must be dried immediately after cleaning and disinfecting.




TBQ GERMICIDAL PRODUCTS ARE NOT RECOMMENDED AS THEY MAY CAUSE DAMAGE TO THE CENTRIFUGE. REFRAIN FROM USING TO PREVENT VOIDING THE WARRANTY.

- Fully/partially halogenated hydrocarbons, ketones, esters, ethers, benzyls, ethyl benzenes, and all other chemicals not prescribed by the manufacturer shall not be used as they may cause damage to the rotor chamber, rotor, tube holders, accessories and centrifuge exterior and void the warranty.

TROUBLESHOOTING

NOTE: The latch must be turned completely clockwise to its stop position for the centrifuge to operate.

<p>The centrifuge does not run</p>	<ul style="list-style-type: none"> ○ Verify that the centrifuge is powered. One of the LED lights should be on. ○ If “Lid not closed” message is displayed, make sure the lid latch is turned completely clockwise to its stop position. ○ If the centrifuge still does not run, contact Customer Service.
<p>The rotor does not spin freely</p>	<ul style="list-style-type: none"> ○ Make sure nothing has fallen into the rotor chamber, following the procedure above. ○ If nothing obstructs the rotor, the rotor may be damaged. Contact Customer Service for further assistance.
<p>The centrifuge makes a rattling noise when running</p>	<ul style="list-style-type: none"> ○ Stop the centrifuge. Open the lid. ○ Wearing PPE, remove tubes and tube holders/buckets and look for fallen objects or debris. Carefully reach inside the rotor chamber with a tool to remove them. ○ Inspect the rotor, tube holders or buckets for damage. ○ If the tube holders or buckets have any damage, even slight, safely dispose of them and replace them. ○ If the rotor appears damaged, contact Customer Service for further assistance.
<p>Excessive noise or vibration when the centrifuge is running</p>	<ul style="list-style-type: none"> ○ Verify that all four centrifuge feet are properly seated on a flat surface. ○ Ensure that the load is balanced according to instructions in the “Balancing Loads” section of this manual. ○ Make sure that nothing has fallen into the rotor chamber.
<p>“Abort” is displayed on the top screen</p>	<p>The centrifugation cycle has been interrupted.</p>
<p>The centrifuge stops and beeps continuously</p>	<p>The load is not balanced. Press the UNLOCK button, open the lid, and balance the load as recommended elsewhere in this manual.</p>
<p>The centrifuge is stuck on one of the settings</p>	<p>Cycle selection is locked. Press the UNLOCK button for 5 seconds.</p>
<p>Only a few cycles can be accessed</p>	<p>The Preset Lock is active. To deactivate it, press the UNLOCK button for 5 seconds, until you hear 2 beeps, then again until the next 2 beeps. All cycles can now be accessed and/or amended.</p>
<p>The cycle time and speed are not set to the desired value</p>	<p>Check the setting by following the instructions in the section on Changing Cycle Settings. If the preset is not the desired length, follow the procedure on the same page to change the preset time.</p>
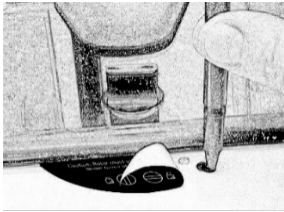
<p>Cycle parameters cannot be changed</p>	<ul style="list-style-type: none"> ○ If cycle selection is locked on one cycle, press the UNLOCK button for 5 seconds. Then, press the GEAR button and follow the instructions elsewhere in this manual. ○ If different saved cycles can be selected but not modified, the centrifuge is in Preset Lock mode. Press the UNLOCK button for 5 seconds until two beeps are heard, then again until the next two beeps. You should now be able to change cycle parameters.
<p>The centrifuge does not unlock after a run is completed</p>	<ul style="list-style-type: none"> ○ Wait until the rotor has come to a complete stop. If the lid knob still cannot be rotated, press the UNLOCK button and try again. <ul style="list-style-type: none"> ○ If no LED light is on, the unit is not powered and the lid will not unlock by conventional means. Remove the latch label and use a pen to manually disengage the locking mechanism. Pull the mechanism towards the control panel and then unlatch and open the lid. ○ If the unit is damaged, contact Customer Service for assistance. 
<p>The lid does not open</p>	<ul style="list-style-type: none"> ○ Ensure that the lid knob is turned fully counterclockwise. ○ If the knob cannot be turned counterclockwise, turn it fully clockwise, press UNLOCK, and turn counterclockwise. ○ If the lid remains locked after this and will not unlock, the electronics may have been damaged. Contact customer service for assistance.
<p>Clicking noise during braking gets loud</p>	<ul style="list-style-type: none"> ○ Make sure that the screw in the center of the rotor is tight.
<p>Lid does not stay up</p>	<ul style="list-style-type: none"> ○ Tighten the center screw on the lid hinge.

GENERAL SPECIFICATIONS

The rotor and accessories are rated for the maximum rotation frequency shown below.

Tube Capacity	24 tubes – 3 to 10 mL 12 tubes – 3 to 15 mL 6 tubes – 50 mL
Dimensions (H x W x D)	9 in x 15 in x 17 in (23 cm x 38 cm x 43 cm)
Weight	31 lbs. (14 kg)
Sound Level	64 dB A
Environmental Range	16 – 32 °C
Voltage	100 -240 VAC
Frequency	50/60 Hz
Power Requirement	415 Watts
Centrifuge Motor	½ H.P. Brushless DC
Max g-Force	3,450 xg
Max Speed	4,500 RPM
Cycle Time	30 sec to 99 min 59 sec (+/- 2%)

CALCULATING THE G-FORCE

The I.F.U.s of tube manufacturers recommend cycles at a minimum G-Force, which can be calculated if you know the RPM and the radius. Use the formula below or go to www.druckerdiagnostics.com/g-force-calculator/.

$$\begin{aligned} &\text{In Centimeters:} \\ \text{RCF or G-force} &= 0.00001118 \times \\ &\text{Rotor Radius (cm)} \times (\text{RPM})^2 \end{aligned}$$

$$\begin{aligned} &\text{In Inches:} \\ \text{RCF or G-force} &= 0.0000284 \times \\ &\text{Rotor Radius (in)} \times (\text{RPM})^2 \end{aligned}$$

Radius 6.0 in (15.3 cm)

REPLACEMENT PARTS

Part No.	Description
7724052K	Foot, rubber (Pack of 4)
02-002-1-0044	Lid Tray Assembly
02-001-0-0008	24-Place Rotor, Horizontal
02-005-1-0012	Motor Assembly
02-006-0-0021	PC Board
7760006	Power cord
03-1-0005-0193	Power Supply, 350W
02-002-1-0037	Lid Assembly
7724071K	Hinge, friction (Pack of 2)
02-002-1-0058	Seal, lid gasket
03-0-0003-0313	Open/Close Label
02-004-0-0012K	4pl Carrier Assembly, 30% CF (Pack of 6)
03-1-0005-0067K	Kit, Lid LED Assembly, Green
03-0-0003-0601	Dash Flex 24 Front Panel Label
00-100-100-010	D/H 12/24 Soft Button & Spacer Replacement Kit
00-100-100-005	Replacement Grommet & Bushing Kit

This operator's manual is part number 03-0-0002-0206 Rev. D

To access previous manuals, please click the relevant link below:

[Revision B](#)

[Revision C](#)

Product Family: DASH Flex Series

Complies with UL61010-1/CSA C22.2 No. 61010-1 and IEC61010-2-020

Protected by U.S. Patents #6,811,531, # 7,422,554, #D718,463, & #D734,489. Other Patents Pending

FDA LISTED



INSTRUCTIONS FOR DISPOSAL OF WEEE BY USERS IN THE EUROPEAN UNION



This product must not be disposed of with other waste. Instead, it is the user's responsibility to dispose of their waste equipment by handing it over to a designated collection point for the recycling of waste electrical and electronic equipment. The separate collection and recycling of your waste equipment at the time of disposal will help to conserve natural resources and ensure that it is recycled in a manner that protects human health and the environment. For more information about where you can drop off your waste equipment for recycling, please contact your local city office, waste disposal service, or where you purchased the product.

Designed, built, and supported in the USA



200 SHADY LANE, SUITE 170 – PHILIPSBURG, PA 16866, USA

+1-866-265-1486 (U.S. ONLY) - +1-814-692-7661

CUSTOMERSERVICE@DRUCKERDIAGNOSTICS.COM

DRUCKERDIAGNOSTICS.COM

