

BOOST

4+ Flex











Operator's Manual



TABLE OF CONTENTS

Symbols.....	3
Caution and Warning Statements.....	4
Model Description	6
Features	6
Intended Use	6
Warranty.....	6
Initial Setup	7
Operation	7
Quick Start	8
Settings	9
Display Cycle Count	9
Review Acceleration Setting	9
Loading.....	10
Carrier Alignment	10
Balancing Loads	11
Care and Preventative Maintenance	11
Cleaning and Disinfection	12
Troubleshooting.....	13
Calculating the G-Force.....	14
General Specifications	15

Symbols

Symbol	Definition	Use
	Caution	Caution to safety hazard. Potential risk of personal injury or damage to the instrument if improperly handled. Consult the manual before proceeding.
	Electrical Hazard	Hazardous Voltage. Potential risk of personal injury or damage to the instrument.
	Manufacturer	Manufacturer of record.
	Electrical and electronic products recycling symbol	Recycle only as electronic waste. Do not dispose in normal waste.
	RoHS Compliant	Compliance with RoHS environmental standards.
	CE Mark	Denotes conformity to specific European directives and regulations.
	UKCA Mark	Denotes conformity to specific UK directives and regulations
	MET Listing	Denotes conformity to specific safety standards and regulations.
	FDA Listed	Denotes that the product has been properly listed with the FDA.
	ISO Certification	Denotes conformity to quality standards and quality management systems.

CAUTION AND WARNING STATEMENTS



This device is intended to be operated by properly trained personnel who have carefully read the operating manual and are familiar with the function of the device. [Refer to the clinical laboratory method specified by the specimen receptacle manufacturer or established by the medical technology for the products applications.]



WARNING: For the safety of both the operator and service personnel, care should be taken when using this centrifuge if handling substances that are known to be toxic, radioactive, or contaminated with pathogenic microorganisms. Use appropriate personal protection equipment (PPE). When Risk Group II materials are used, (as identified in the World Health Organization “Laboratory Bio-Safety Manual”), a Bio-Seal should be employed. If materials of a higher risk group are being used, more than one level of protection must be provided. The use of flammable or explosive materials as well as those materials which have a vigorous chemical reaction is prohibited.



Users of centrifuge should validate the processing of their disposable for their specific application prior to use.



Operation of this equipment in a manner not specified by the manufacturer may impair the protection provided by the equipment.



Unplug the centrifuge before cleaning or performing maintenance.



WARNING: Do not make modifications to or remove any hardware from rotor without prior authorization from Drucker Diagnostics.



WARNING: Only use Drucker Diagnostics components in this centrifuge



WARNING: Inspect centrifuge for cracks or physical damage to cabinet, lid, rotor, or tube holders. Damage may result in unsafe operation. Discontinue use until repairs have been performed.



For your safety and the durability of the machine, never transport or store centrifuge with tube holders inside machine.



Electrical safety protection is provided by properly connecting the centrifuge to earth ground. Use only the manufacturer provided line cord and ensure that it is connected to a properly grounded power receptacle. Failure to do so will result in an electrical hazard.



This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with this operator manual, may cause interference to radio communications.



Operation of this equipment in a residential area may cause interference, in which case the user will be required to correct the interference at his own expense.



Due to the lack of the possibility of human exposure, all Drucker centrifuges and accessories sold by Drucker Diagnostics, Inc. are compliant without any special labeling required by the California Safe Drinking Water and Toxic Enforcement Act (Proposition 65).



The maximum combined weight allowed to be loaded into each carrier (including Drucker provided components) shall not exceed maximum specification. Refer to **General Specifications** section for maximum carrier load specification.

MODEL DESCRIPTION

Satisfy your lab's most diverse processing requirements with the Boost 4+ Flex. Customize settings and check every detail on the digital display. Program up to 10 custom cycles.

This general-purpose laboratory centrifuge may also be used to spin approved containers with biologics, chemicals (non-flammable, non-explosive, non-volatile, and non-highly reactive), and environmental samples.

FEATURES

- Preset cycles have been programmed for time, speed, acceleration and braking to facilitate matching validated cycles and Ortho Ascent IFUs.
- A digital cycle counter tracks the number of cycles the centrifuge has run.
- Lid lighting indicates the centrifuge's status (ready, running, done), informing the operator when tubes are ready for the analyzer and preventing tubes from being left in the centrifuge longer than necessary (patent pending).
- A traditional audible alert indicates the completion of the cycle. The audible alert can be muted.
- Cool-Flow design prevents overheating of samples by using ambient air to keep specimens at room temperature.
- The carriers are fiber reinforced for high strength, durability, and years of trouble-free use.
- A clear lid permits safe observation of samples and optical calibration of speed.
- The lid safety system prevents the centrifuge from operating unless the lid is closed and latched.
- The lid safety system only allows entry into the centrifuge after the rotor has completely stopped.
- The high-power brushless motor provides years of operation with no routine maintenance.

INTENDED USE

General purpose laboratory centrifuge, intended for the density-based separation of fluids through centripetal acceleration.

WARRANTY

Drucker Diagnostics warrants that this centrifuge is free from defects in workmanship and parts for 2 years.

INITIAL SETUP



For your safety, only setup the centrifuge within the appropriately designated operating environment. Refer to **General Specifications** section for pertinent environmental condition requirements when using the centrifuge.

- Unpack and verify that all the following are included:
 - Centrifuge
 - Power cord
 - Carriers
 - Ortho Ascent accessories
 - Quick Start Insert or Operator's Manual
- Place the centrifuge on a flat and level surface. A bench top clearance height of 22" (56 cm) is required to open the lid.
- The centrifuge should have 6" (15 cm) of clear space around the centrifuge. Proper ventilation is necessary to prevent the overheating of samples as well as premature failure of the centrifuge. Choose an area which allows unencumbered air flow.
- No hazardous material shall be permitted in the clearance envelope during operation.
- The operator time within the envelope shall be limited to the time necessary for loading, unloading, and centrifuge operation only.
- Plug the line cord into the centrifuge.
- Plug the line cord into an approved electrical outlet.



Be sure that the electrical outlet is always accessible, as the line cord is the means of emergency disconnection.

OPERATION

- Place disposable(s) into the carriers, utilizing inserts as appropriate.



The maximum combined weight allowed to be loaded into each carrier (including Drucker provided components) shall not exceed maximum specification. Refer to **General Specifications** section for maximum carrier load specification.




Be sure to follow the rules for balanced loads. Refer to **Balancing Loads** section.

- Close the lid and turn the lid knob clockwise to its complete stop position.
- The digital screen shows the currently selected cycle. To select another cycle, press the UP or DOWN button in succession until the desired cycle is selected.
- Pushing the START button on the control panel will start the spin cycle.

- When the cycle is completed, the rotor will slow to a complete stop and the lid light will flash.
- The unlocking mechanism will engage for 60 seconds allowing entry into the rotor chamber. To unlock after more than 60 seconds have elapsed, press the UNLOCK button. The lid will unlock for another 15 seconds.
- Turn the lid knob counterclockwise and open the lid. The lid light will turn off.
- You may now safely remove the disposable(s).

QUICK START


The top screen display alternates between the name of the currently selected cycle and speed. The bottom screen displays the time setting.

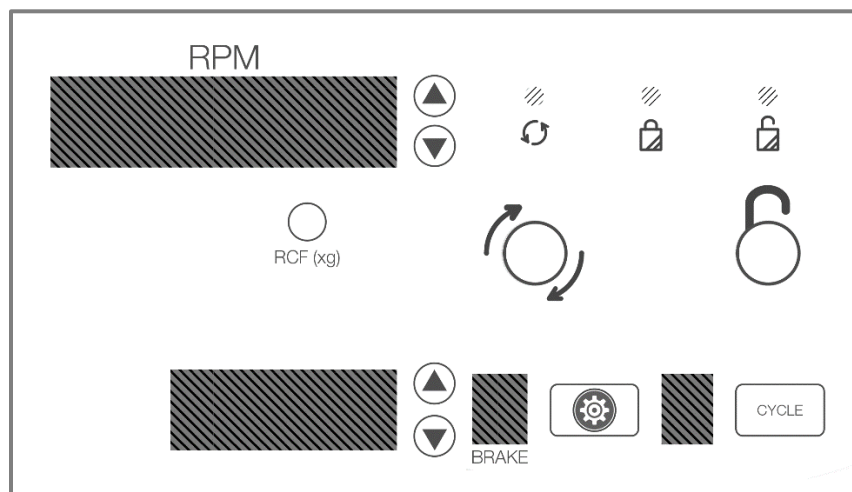


	Start	Begins running the cycle displayed on the screen. The lid must be closed.
	Unlock	Allows access into the rotor chamber by engaging the unlocking mechanism. Entry is only possible when the rotor is stopped.
	Stop	Pressing the UNLOCK button during operation will terminate the run and unlock the lid after the rotor has come to a stop.
	Cycle Selection	Press the CYCLE button to select the desired saved cycle.

SETTINGS

DISPLAY CYCLE COUNT

 [HOLD]	Display Cycle Count	With the lid open and the unit powered, press and hold the START button for five (5) seconds. The cycle count will be displayed until the START button is released.
--	----------------------------	---

REVIEW ACCELERATION SETTING

 [HOLD]	Display Cycle Count	With the lid open and the unit powered, press and hold the GEAR button. The TIME display will show "ACCL", and the BRAKE display will show the pertinent acceleration setting, until the GEAR button is released.
--	----------------------------	--



LOADING

CARRIER ALIGNMENT

To ensure proper operation, align carriers and rotor utilizing the arrows indicated below.

Verify that carriers are hanging from the rotor pins and swing freely into the horizontal position.

Rotor P/N 03-1-0001-0138

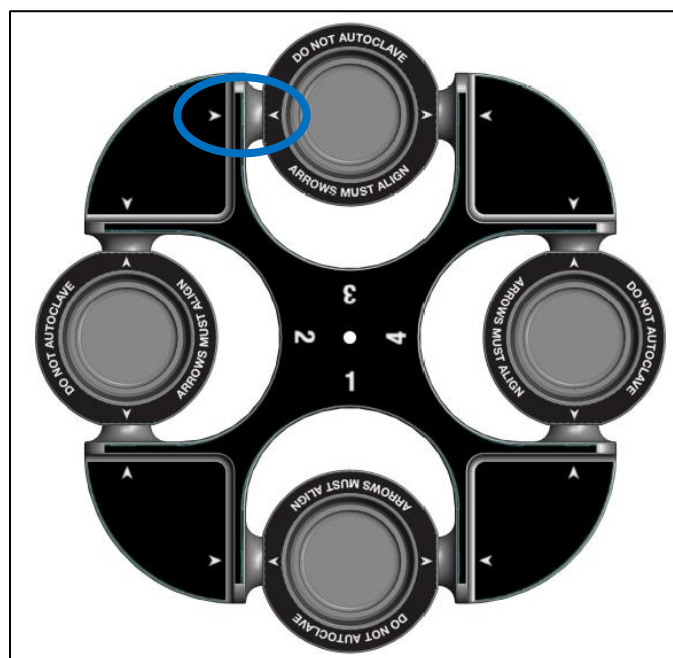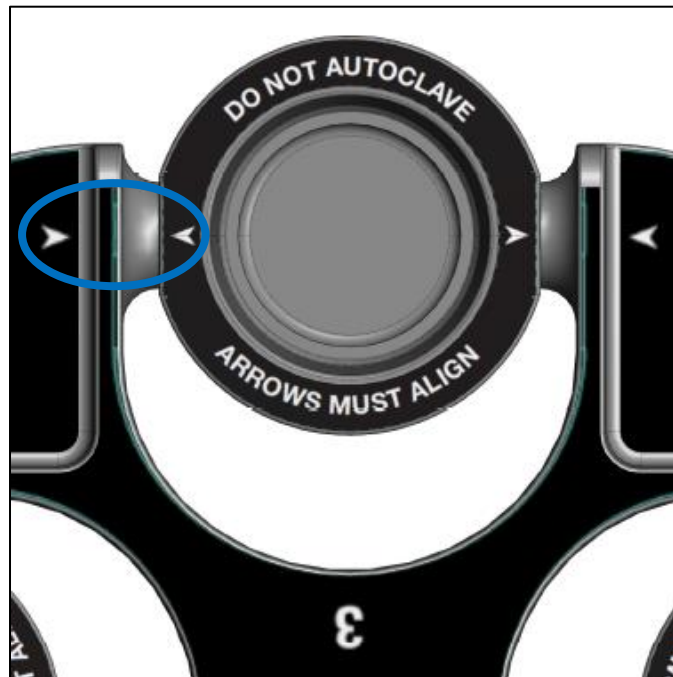

BALANCING LOADS


Your centrifuge must contain a balanced load to work properly. Spinning balanced loads will extend the life of the centrifuge and produce better results.



The amount of imbalance allowed to be loaded into each opposing carrier shall not exceed maximum specification. Refer to **General Specifications** section for maximum imbalance specification.

Use the following rules when loading the rotor. If an odd number of samples is to be spun, fill a tube with water to match the weight of the unpaired sample and place it across from this sample.

Opposing carriers must be equally loaded, empty, or loaded with equally weighted samples.

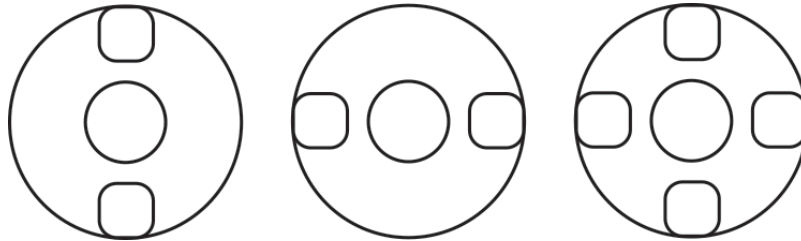

CARE AND PREVENTATIVE MAINTENANCE

With proper care and maintenance, your centrifuge will provide years of laboratory service. For proper care, the following steps should be taken:

- **Always Spin Balanced Loads:** Make certain that you are always spinning a balanced load, as shown in the previous section. These centrifuges have a unique counter balanced motor mounting design which produces excellent vibration dampening. However, out-of-balance loads may break glass test tubes and may produce unsatisfactory separation results. Proper load balancing will improve sample separation and extend the life of the centrifuge.
- **Motor and Electrical Maintenance:** The highest quality electrical components have been selected for the centrifuges and should not need maintenance or servicing for the life of the centrifuge.
- **Carrier Replacement:** It is recommended that the carriers be replaced after 24 months of use. Inspect carriers regularly for cracks. If cracks are discovered, replace immediately.
- **Remove Accessories Before Moving:** All carriers, inserts, and disposables must be removed from the rotor chamber before transporting or storing the centrifuge to prevent damage and injury.

CLEANING AND DISINFECTION

To prolong the life of the centrifuge, cleaning and disinfection is recommended every six months or whenever there is a spillage or tube breakage. Contaminants must be removed immediately, or corrosion and premature degradation of components can occur. Before using any cleaning or decontamination methods other than those recommended by the manufacturer, users should verify with the manufacturer that the proposed method will not damage the equipment.

- Unplug the centrifuge before cleaning.
- Use appropriate personal protective equipment (PPE).
- Apply cleaning solutions with a dampened towel or cloth ONLY. Do not spray or pour cleaning solution directly onto or into the centrifuge. Do not saturate the centrifuge or submerge the centrifuge in water or other cleaning solutions as this will cause damage, create a safety risk, and void the warranty.
- ONLY isopropyl alcohol or a 10% (5500 PPM) bleach solution should be used to disinfect the centrifuge and its accessories.
- All surfaces must be dried immediately after cleaning and disinfecting.
- Fully/partially halogenated hydrocarbons, ketones, esters, ethers, benzyls, ethyl benzenes, and all other chemicals not prescribed by the manufacturer shall not be used as they may cause damage to the rotor chamber, rotor, carriers, accessories, and centrifuge exterior and void the warranty.




TBQ Germicidal products are not recommended as they may cause damage to the centrifuge. Refrain from using to prevent voiding the warranty.

TROUBLESHOOTING

NOTE: The latch must be turned completely clockwise to its stop position for the centrifuge to operate.

<p>The centrifuge does not run</p>	<ul style="list-style-type: none"> ○ Verify that the centrifuge is powered. One of the LED lights should be on. ○ If “LID” message is displayed, make sure the lid latch is turned completely clockwise to its stop position. ○ If the centrifuge still does not run, contact Customer Service.
<p>The rotor does not spin freely</p>	<ul style="list-style-type: none"> ○ Make sure nothing has fallen into the rotor chamber, following the procedure above. ○ If nothing obstructs the rotor, the rotor may be damaged. Contact Customer Service for further assistance.
<p>The centrifuge makes a rattling noise when running</p>	<ul style="list-style-type: none"> ○ Stop the centrifuge. Open the lid. ○ Wearing PPE, remove tubes and carriers/inserts and look for fallen objects or debris. Carefully reach inside the rotor chamber with a tool to remove them. ○ Inspect the rotor, carriers or inserts for damage. ○ If the carriers or inserts have any damage, even slight, safely dispose of them and replace them. ○ If the rotor appears damaged, contact Customer Service for further assistance.
<p>Excessive noise or vibration when the centrifuge is running</p>	<ul style="list-style-type: none"> ○ Verify that all four centrifuge feet are properly seated on a flat surface. ○ Ensure that the load is balanced according to instructions in the “Balancing Loads” section of this manual. ○ Make sure that nothing has fallen into the rotor chamber.
<p>“Abort” is displayed on the top screen</p>	<p>The centrifugation cycle has been interrupted.</p>
<p>The centrifuge stops and beeps continuously</p>	<p>The load is not balanced. Press the UNLOCK button, open the lid, and balance the load as recommended elsewhere in this manual.</p>
<p>The centrifuge is stuck on one of the settings</p>	<p>Cycle selection is locked. Press the UNLOCK button for 5 seconds.</p>
<p>Only a few cycles can be accessed</p>	<p>The Preset Lock is active. To deactivate it, press the UNLOCK button for 5 seconds, until you hear 2 beeps, then again until the next 2 beeps. All cycles can now be accessed and/or amended.</p>
<p>The cycle time and speed are not set to the desired value</p>	<p>Check the setting by following the instructions in the section on Setting or Modifying a Saved Cycle. If the preset is not the desired length, follow the procedure on the same page to change the preset time.</p>
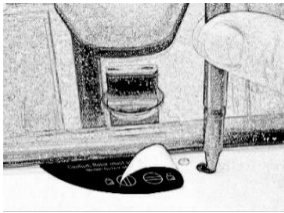
<p>Cycle parameters cannot be changed</p>	<ul style="list-style-type: none"> ○ If cycle selection is locked on one cycle, press the UNLOCK button for 5 seconds. Then, press the GEAR button and follow the instructions elsewhere in this manual. ○ If different saved cycles can be selected but not modified, the centrifuge is in Preset Lock mode. Press the UNLOCK button for 5 seconds until two beeps are heard, then again until the next two beeps. You should now be able to change cycle parameters.
<p>The centrifuge does not unlock after a run is completed</p>	<ul style="list-style-type: none"> ○ Wait until the rotor has come to a complete stop. If the lid knob still cannot be rotated, press the UNLOCK button and try again. <ul style="list-style-type: none"> ○ If no LED light is on, the unit is not powered and the lid will not unlock by conventional means. Remove the latch label and use a pen to manually disengage the locking mechanism. Pull the mechanism towards the control panel and then unlatch and open the lid. ○ If the unit is damaged, contact Customer Service for assistance. 
<p>The lid does not open</p>	<ul style="list-style-type: none"> ○ Ensure that the lid knob is turned fully counterclockwise. ○ If the knob cannot be turned counterclockwise, turn it fully clockwise, press UNLOCK, and turn counterclockwise. ○ If the lid remains locked after this and will not unlock, the electronics may have been damaged. Contact customer service for assistance.
<p>Lid does not stay up</p>	<ul style="list-style-type: none"> ○ Tighten the center screw on the lid hinge.

CALCULATING THE G-FORCE

The I.F.U.s of tube manufacturers recommend cycles at a minimum G-Force, which can be calculated if you know the RPM and the radius. Use the formula below or go to www.druckerdiagnostics.com/g-force-calculator/.

In Centimeters:
 $RCF \text{ or } G\text{-force} = 0.00001118 \times \text{Rotor Radius (cm)} \times (\text{RPM})^2$

In Inches:
 $RCF \text{ or } G\text{-force} = 0.0000284 \times \text{Rotor Radius (in)} \times (\text{RPM})^2$

Radius 6.46 in (16.4 cm)
(Displayed RCF is calculated at 16.4 cm radius,
and is rounded to the nearest 25 xg.)

GENERAL SPECIFICATIONS

The rotor and accessories are rated for the maximum rotation frequency shown in the table below.

Disposable Capacity	4 disposables (Up to diameter of 4.6 cm / 1.81 in) (Up to length of 13.0 cm / 5.12 in)
Dimensions (H x W x D)	15 in x 17 in x 10 in (38 cm x 43 cm x 25 cm)
Weight	39 lbs. (17 kg)
Sound Level	75 dB A
Supply Voltage (+/- 10%)	100-240 VAC
Supply Frequency	50-60 Hz
Current Consumption	4.2 amps at 115 VAC 2.1 amps at 230 VAC
Centrifuge Motor	½ H.P. Brushless
Max Speed	4,000 RPM (+/- 100)
Max G-Force	2,950 xg
Max Carrier Load	140.0 grams
Cycle Time	.5 to 99 minutes (+/- 2%)
Environmental Conditions	
Set-up Site	Indoor Use Only
Altitude	Up to 2,000m from Sea Level
Ambient Temperature	5 °C to 40 °C
Humidity	Maximum relative humidity 80% for temperatures up to 31 °C, decreasing linearly to 50% relative humidity at 40 °C
Overtoltage Category	II
Pollution Degree	2

This operator's manual is part number 03-0-0002-0307, Rev. B

To access previous manuals, please click relevant link below.

[Revision A](#)

Product Family: BOOST Series (BOOST 4+ Flex)

Complies with UL61010-1/CSA C22.2 No. 61010-1 and IEC61010-2-020

Protected by U.S. Patents #6,811,531, # 7,422,554, #D718,463, & #D734,489. Other Patents Pending

FDA LISTED



INSTRUCTIONS FOR DISPOSAL OF WEEE BY USERS IN THE EUROPEAN UNION



This product must not be disposed of with other waste. Instead, it is the user's responsibility to dispose of their waste equipment by handing it over to a designated collection point for the recycling of waste electrical and electronic equipment. The separate collection and recycling of your waste equipment at the time of disposal will help to conserve natural resources and ensure that it is recycled in a manner that protects human health and the environment. For more information about where you can drop off your waste equipment for recycling, please contact your local city office, waste disposal service, or where you purchased the product.

Designed, built, and supported in the USA



Drucker Diagnostics
200 SHADY LANE, SUITE 170 – PHILIPSBURG, PA 16866, USA
+1-877-231-3115 (U.S. ONLY) - +1-814-692-7661
CUSTOMERSERVICE@DRUCKERDIAGNOSTICS.COM
DRUCKERDIAGNOSTICS.COM

