



614VET VETSCAN IMAGYST™ READY CENTRIFUGE

Process Vetscan Imagyst™ specimen tubes at the push of a button. The 614VET makes it easy to centrifuge at any veterinary practice.

Scan to buy or learn more:



It's designed for the unique needs of you and your team - but it's built on the technology of the world's best selling reference lab centrifuge. That means years of engineering to deliver maximum durability, reliability, and value for your practice.

The Easiest Way to Prepare Imagyst™ Samples

No guesswork - just insert your prepared Imagyst™ tubes and start the cycle. The 614VET is optimized for Imagyst™ with a single preset runtime of 2 minutes at 370 xg.

Conserve Your Counter Space

It's easy to find space for the 614VET. This compact unit only requires 1 square foot of counter. And it's easy to move and store if needed when not in use.

Unmatched Reliability

The 614VET is a workhorse. It will run for years with no maintenance and no downtime. It's backed by a two year unlimited warranty and a lifetime warranty on the rotor.

A Great Value

The 614VET fits your budget. The unit arrives ready to run Imagyst™ sample tubes with rotor, tube holders, and 2-year warranty included in the affordable price.

Capacity	Applications	Centrifugation	Included	Warranty
6 Imagyst™ tubes or 15 mL/125 mm	Imagyst™ tubes General veterinary	Fixed Angle 370 xg Max RCF 1,750 Max RPM	Rotor Tube Holders	2 years Rotor: lifetime warranty Extended warranty available

Designed, built, and supported in the USA

614VET Centrifuge *Preset for Imagyst™*


Scan to buy:



Safety lid automatically locks during operation and will not open until the centrifuge has stopped

Reinforced guard bowl for maximum safety

Brushless motor ensures quiet, maintenance-free operation

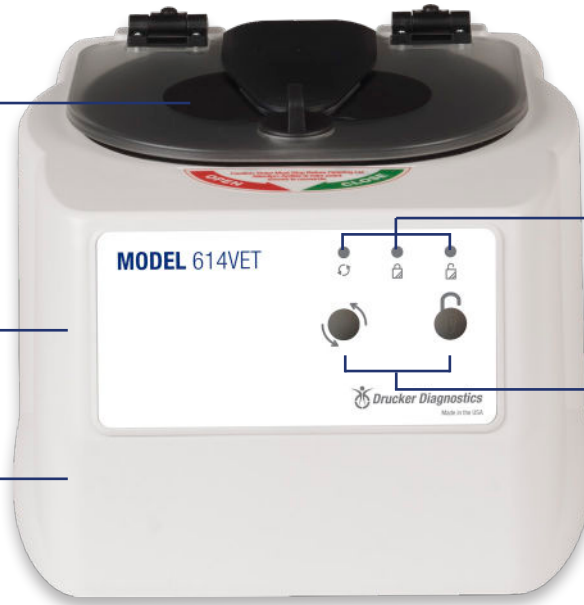

Lights illuminate to indicate running, locked, unlocked

Simple two-button controls

Part Numbers	
Centrifuge	00-118-226-000
75-100 mm tube holder	7713079 (1 qty)
125 mm tube holder	7713032 (1 qty)
SmartView Platform	00-079-009-001
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Specifications	
Sound level	64 dBA
Dimensions	11.75 x 14 x 8.75 in 30 x 35.5 x 22 cm
Weight	11 lbs (5 kg)
Max Acceleration	10 seconds
Max Deceleration	20 seconds
Power	115 V – 50/60 Hz

FDA LISTED



RoHS
Compliant



ISO13485
certified



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661

