


642E **SINGLE-SPEED ELECTRONIC CENTRIFUGE**

Part # 00-018-009-001

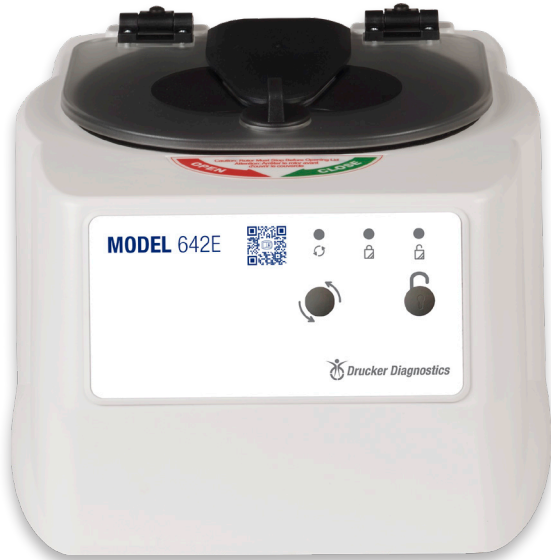

The 642E is the **world's best-selling clinical laboratory centrifuge**. Easy two-button operation ensures repeatability and standardization, while advanced design delivers high-quality blood separation every time. Its preset speed and runtime make the 642E the ideal choice for consistency in doctor's offices, patient service centers, and outreach locations.

Easy to Use

Spin samples with the press of just one button. The 642E processes tubes up to 16 x 100 mm or 10 mL for chemistry, coag, and PPP.

Simple Settings

The 642E features a preset runtime of 10 minutes at 1600 xg (3380 RPM) for all common chemistry needs.

Built to Last

No routine maintenance. The brushless motor and unique convection-based cooling system ensure years of high-quality separations.

Engineered for Safety

Protect your lab with the clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Get the world's #1 clinical laboratory centrifuge at a price that fits any budget. The unit arrives ready for use with the horizontal rotor, carriers, and a 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

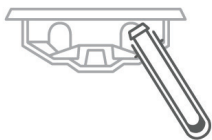


Rotation Angle

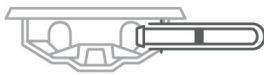

Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag, PPP	3380 Max RPM 1600 xg Max RCF	6 x up to 16 x 100 mm or 10 mL tubes	Horizontal Rotor 6 x 16 x 100 mm Blue Carrier	2 years Extended warranty options available

Designed, built, and supported in the USA

642E Centrifuge

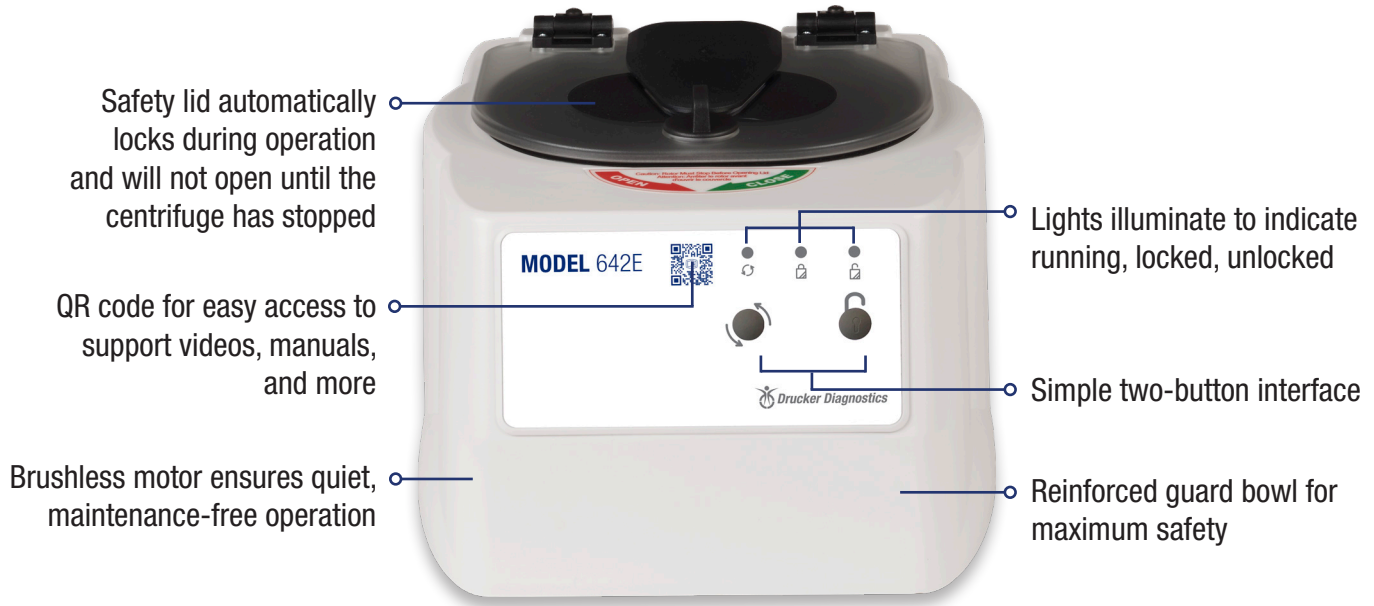

Standard Preset Cycles	Time	RPM	G-Force
Cycle 1	10 min	3380	1600

Part Numbers

	Centrifuge	00-018-009-001
Included	16 x 100 mm Carrier, Blue (6 qty)	7713079K
	1-1.5 mL Microtainer Tube Adapter (6 qty)	03-1-0007-0059K
Optional	1 yr Warranty Extension	03-0-0007-0001
	2 yr Warranty Extension	03-0-0007-0002
	3 yr Warranty Extension	03-0-0007-0006

(adapters available for a variety of tube sizes)

Specifications

Sound level	64 dBA
Dimensions (W x D x H)	12 x 14 x 9 in 30.5 x 35.6 x 22.9 cm
Weight	11 lbs (5 kg)
Max RPM	3380 RPM
Max G-Force	1600 xg
Power	115 V – 50/60 Hz



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661