

BOOST 2+ FLEX **ADJUSTABLE CAPACITY CENTRIFUGE**

Part # 00-783-009-000

Spin large tubes in a small centrifuge with the BOOST 2+ Flex. Fitted with our largest buckets to hold tubes up to 1.8"W x 5.1"H (4.6 x 13.0 cm). A lightweight fit for any space with a small footprint of 13 x 15 x 9 inches.

Fully Customizable

The BOOST 2+ Flex offers a variety of extensive customization options. Personalize the front panel label and LED lid lights to match your unique branding, or pre-program custom RPM and g-force settings to meet your preferred processing specifications.

Spins Your Cycles

Program and recall up to 10 preset cycles with the BOOST 2+ Flex's backlit digital display. Can be set to run preset cycles only, or may be locked to just one preset cycle if desired.

Fits Your Tubes

Experience unmatched convenience with BOOST 2+ Flex's unique insert system. Simply select the correct adapter for your tube sizes and switch them out as needed to spin every tube.

Safe and Durable

Enjoy years of worry-free operation with BOOST's imbalance protection and shatterproof lid, which can only be opened when the rotor is stopped - even if the unit loses power. Rigorous electrical testing ensures our centrifuges meet hospital equipment safety standards.

A Great Value

Unbeatable quality and flexibility for any budget. Each BOOST 2+ Flex centrifuge is ready to run with swing-out horizontal rotor, buckets, and 2-year warranty included in the affordable price.




A Perfect Fit For Any Tube

Our patented rotor and suite of custom-designed adapters ensure a perfect fit for any tube or disposable.

Applications	Max Speed	Max Capacity	Included	Warranty
Variable Tube Sizes Custom Disposables	4400 RPM 3000 xg RCF	2 x Up to 1.9 x 5.4 in (4.8 x 13.7 cm) or 140 grams	Horizontal Rotor 2 x BOOST Bucket for Molded Rotor	2 years Extended warranty options available

BOOST 2+ Flex

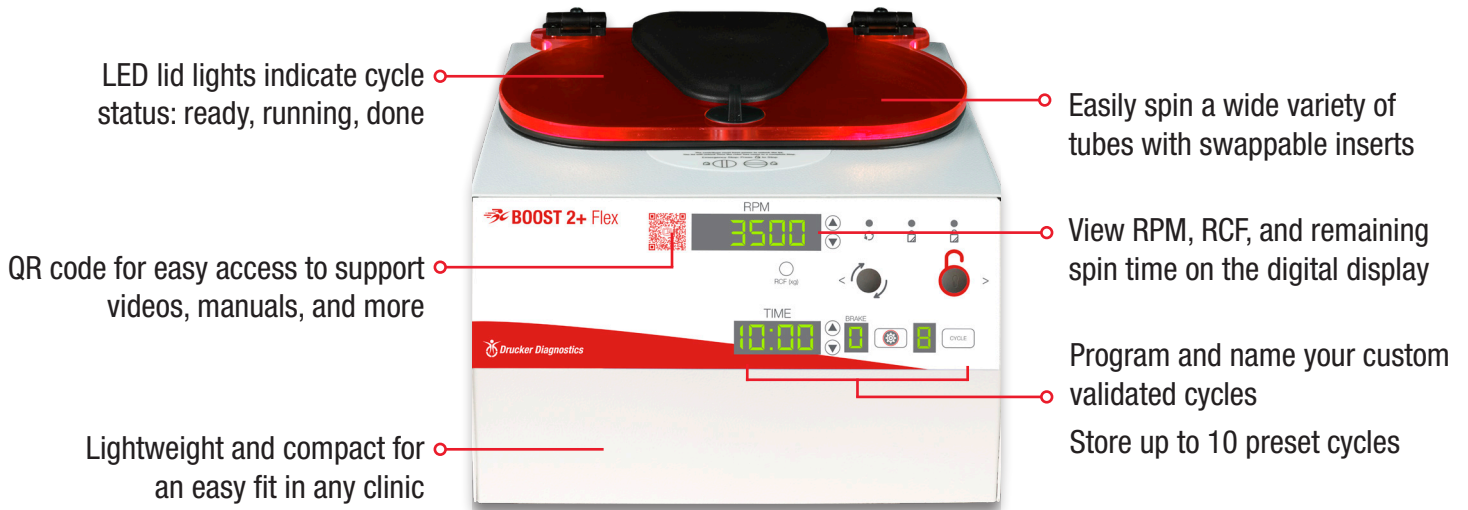

Free 30-Day Evaluation (U.S. Only)

Part Numbers		Specifications		
	Centrifuge	00-783-009-000	Sound Level	70 dBA
Included	BOOST Bucket for Molded Rotor (2 qty)	03-1-0007-0123K	Max Speed RPM	4400 RPM
	.68" BOOST Bucket Adapter, Red (2 qty)*	03-1-0007-0112K	Max G-Force	3000 xg
Optional	.85" BOOST Bucket Adapter, Red (2 qty)*	03-1-0007-0113K	Dimensions (W x D x H)	13 x 15 x 9 in 33 x 38.1 x 22.9 cm
	1.03" BOOST Bucket Adapter, Red (2 qty)*	03-1-0007-0114K	Weight	34 lbs (15.42 kg)
	1.2" BOOST Bucket Adapter, Red (2 qty)*	03-1-0007-0115K	Power	100 - 240V - 50/60 Hz
	1.34" BOOST Bucket Adapter, Red (2 qty)*	03-1-0007-0116K	Resting Tube Angle	20° (from vertical)
	1.56" BOOST Bucket Adapter, Red (2 qty)*	03-1-0007-0117K	Tube Spin Angle	90° (horizontal)
	1 yr Warranty Extension	03-0-0007-0003	Maximum Load	140 g per bucket
	2 yr Warranty Extension	03-0-0007-0004		
3 yr Warranty Extension	03-0-0007-0005			

*also available in black adapters



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661