

DASH APEX 12 SET & LOCK STAT CENTRIFUGE

Part # 00-083-009-003

Discover a better way to spin STAT with Drucker's DASH Apex 12 STAT centrifuges. Minimize lab error, reduce turnaround time, and streamline training with the DASH Apex's simple controls and powerful performance.



Easy to Monitor

Prevent forgotten samples and reduce TAT with LED lid lighting. Bright lights indicate when the centrifuge is ready to load (off), running (on), or done (flashing). Lights continue flashing until the lid is opened.

Easy to Use

The simple two-button interface makes it easy to select a preset or program your custom validated STAT cycles. Use all three presets, or set and lock to a single cycle for maximum repeatability. Spin tubes up to 16 x 100 mm or 10 mL with the same rotor and carriers.

Eliminate Wait Time

Further shorten TAT by using a centrifuge array in accessioning to eliminate queuing. The array takes advantage of the DASH Apex 12's compact size to allow samples to be spun as soon as they arrive.

Safe and Durable

Enjoy years of worry-free operation with the DASH Apex 12's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each centrifuge comes ready to run with a horizontal rotor, carriers, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

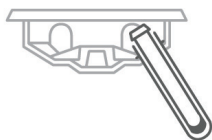

Rotation Angle

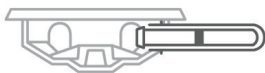

<p>Applications</p> <p>STAT Chemistry, Coag, PPP</p>	<p>Max Speed</p> <p>5200 Max RPM 4000 xg Max RCF</p>	<p>Max Capacity</p> <p>12 x up to 16 x 100 mm or 10 mL</p>	<p>Included</p> <p>Horizontal Rotor 12 x 16 x 100 mm CF Carrier</p>	<p>Warranty</p> <p>2 years Extended warranty options available</p>
---	---	---	--	---

DASH APEX 12



Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles	Time	RPM	G-Force
Cycle 1	3 min	5200	4000
Cycle 2	5 min	5200	4000
Cycle 3	7 min	3600	2000

Part Numbers	
Centrifuge	00-083-009-003
Included	16 x 100 mm Carrier, Carbon Fiber (12 qty) 03-1-0007-0046
	1-1.5 mL Microtainer Tube Adapters (6 qty) 03-1-0007-0059K
Optional	1 yr Warranty Extension 03-0-0007-0003
	2 yr Warranty Extension 03-0-0007-0004
	3 yr Warranty Extension 03-0-0007-0005

Specifications	
Sound Level	61 dBA
Max RPM	5200 RPM
Max G-Force	4000 xg
Dimensions (W x D x H)	13 x 15 x 9 in 33 x 38.1 x 22.9 cm
Weight	30 lbs (13.7 kg)
Power	100 - 240 (+/- 10%) - 50/60 Hz

(adapters available for a variety of tube sizes)



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661