

## DASH APEX 24 SET & LOCK STAT CENTRIFUGE

Part # 00-084-009-005

Discover a better way to spin STAT with Drucker's DASH Apex 24 STAT centrifuges. Minimize lab error, reduce turnaround time, and streamline training with the DASH Apex's simple controls and powerful performance.

### Easy to Monitor

Prevent forgotten samples and reduce TAT with LED lid lighting. Bright lights indicate when the centrifuge is ready to load (off), running (on), or done (flashing). Lights continue flashing until the lid is opened.

### Easy to Use

The simple two-button interface makes it easy to select a preset or program your custom validated STAT cycles. Use all three presets, or set and lock to a single cycle for maximum repeatability.

### High Capacity STAT

The DASH Apex 24 can spin 24 tubes up to 16 x 100 mm or 10 mL, 12 tubes up to 16 x 125 mm or 15 mL, and 6 tubes up to 50 mL using a variety of different carriers and adapters.

### Safe and Durable

Enjoy years of worry-free operation with the DASH Apex 24's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.

### A Great Value


Each centrifuge comes ready to run with a horizontal rotor, carriers, and 2-year warranty included in the affordable price.




### 2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

#### Loading Angle

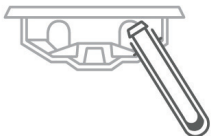

#### Rotation Angle

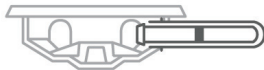

Applications	Max Speed	Max Capacity	Included	Warranty
STAT Chemistry, Coag, PPP	4300 Max RPM 3150 xg Max RCF	24 x up to 16 x 100 mm or 10 mL  12 x up to 16 x 125 mm or 15 mL  6 x 50 mL	Horizontal Rotor  6 x 16 x 125 mm 4-Place CF Carrier	2 years  Extended warranty options available

**Designed, built, and supported in the USA**

# DASH APEX 24



## Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles	Time	RPM	G-Force
Cycle 1	5 min	4300	3150
Cycle 2	6 min	4300	3150
Cycle 3	7 min	3500	2100

Part Numbers	
Centrifuge	00-084-009-005
Included	16 x 125 mm 4-Place Carrier, Carbon Fiber (6 qty)
	02-004-0-0012K
Optional	50 mL Carrier, Carbon Fiber (6 qty)
	02-004-0-0014K
	1-1.5 mL Tube Adapter, Purple (6 qty)
	7713065
	1 yr Warranty Extension
03-0-0007-0003	
2 yr Warranty Extension	
03-0-0007-0004	
3 yr Warranty Extension	
03-0-0007-0005	

*(adapters available for a variety of tube sizes)*

Specifications	
Sound Level	64 dBA
Max RPM	4300 RPM
Max G-Force	3150 xg
Dimensions (W x D x H)	15 x 17 x 9 in 38.1 x 43.2 x 22.9 cm
Weight	39 lbs (17 kg)
Power	100 - 240 (+/- 10%) - 50/60 Hz



druckerdiagnostics.com  
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)  
+1-814-692-7661