

DASH APEX 4 COMPACT STAT CENTRIFUGE

Part # 00-053-009-000

The DASH Apex 4 is a high-speed, fixed angle benchtop centrifuge designed to rapidly separate fluids. Its compact size saves on bench space and allows samples to be moved easily around the lab.

Compact Size

Weighing only 4.3 lbs (1.95 kg) and measuring 7 x 7 x 7 in (17.8 x 17.8 x 17.8 cm), the DASH Apex 4 is the smallest model in the DASH STAT centrifuge family.

Easy to Use

The simple two-button interface makes it easy to select a preset. Use all three presets, or set and lock to a single cycle for maximum repeatability.

Reduce Turnaround Time

Further shorten TAT by using a centrifuge array in accessing to eliminate queuing. The array takes advantage of the DASH Apex 4's compact size to allow samples to be spun as soon as they arrive.

Safe and Durable

Enjoy years of worry-free operation with the DASH Apex 4's brushless motor, reinforced carbon fiber components, and a two-year calibration guarantee. Enhance lab safety with imbalance protection and a shatterproof lid, which locks until the centrifuge has stopped completely and opens.

A Great Value


Start spinning samples right out of the box, no assembly or additional accessories required. Each centrifuge comes ready-to-run with a fixed angle rotor, adapters, and 2-year warranty included in the affordable price.




Fixed Angle Rotor

The Apex 4's fixed angle rotor is compatible with recommended IFUs for all common tube manufacturers, including Greiner and BD. Spin four tubes up to 5 mL, or 1.5-2 mL micro BD Microtainers.

Loading Angle

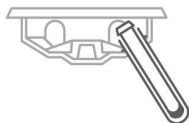

Rotation Angle

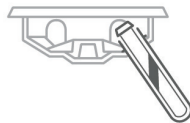

Applications	Max Speed	Max Capacity	Included	Warranty
STAT Chemistry, Coag, PPP	8500 Max RPM 4440 xg Max RCF	4 x 13 x 75 mm or 5 mL	4-Place Fixed Angle Rotor 4 x 13 x 75 mm Tube Adapter	2 years Extended warranty options available

Designed, built, and supported in the USA

DASH APEX 4



Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles	Time	RPM	G-Force
Cycle 1	30 sec	8500	4440
Cycle 2	2 min	8500	4440
Cycle 3	3 min	8500	4440

Part Numbers	
Centrifuge	00-053-009-000
Included	
13 x 75 mm Adapter (4 qty)	7713064-4K
1-1.5 mL Microtainer Tube Adapters (6 qty)	03-1-0007-0059K
Optional	
1 yr Warranty Extension	03-0-0007-0003
2 yr Warranty Extension	03-0-0007-0004
3 yr Warranty Extension	03-0-0007-0005

(adapters available for a variety of tube sizes)

Specifications	
Sound Level	59 dBA
Max RPM	8500 RPM
Max G-Force	4440 xg
Dimensions (W x D x H)	7 x 7 x 7 in 17.8 x 17.8 x 17.8 cm
Weight	4.3 lbs (1.95 kg)
Power	100 – 240 (+/- 10%) VAC – 50/60 Hz



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661