

## DASH APEX 4 SET-AND-LOCK STAT CENTRIFUGE

Part # 00-053-009-000

The DASH Apex 4 is a high-speed benchtop centrifuge designed to rapidly separate fluids in a syringe-type collection tube. Its compact size saves on bench space and allows samples to be moved easily around the lab.

### Compact Size

Weighing only 4.3 lbs (1.95 kg) and measuring 7 in x 7 in x 7 in (18 cm x 18 cm x 18 cm), the DASH Apex 4 is the smallest model in the DASH STAT centrifuge family.

### Easy to Use

The simple two-button interface makes it easy to select a preset. Use all three presets, or set and lock to a single cycle for maximum repeatability.

### Reduce Turnaround Time

Eliminate wait time by using a centrifuge array in accessioning to eliminate queuing. The array takes advantage of the DASH Apex 4's compact size to allow samples to be spun as soon as they arrive.

### Safe and Durable

Enjoy years of worry-free operation with the DASH Apex 4's brushless motor, reinforced carbon fiber components, and a two-year calibration guarantee. Enhance lab safety with imbalance protection and a shatterproof lid, which locks until the centrifuge has stopped completely.

### A Great Value


Start spinning samples right out of the box, no assembly or additional accessories required. Each centrifuge comes ready-to-run with a fixed angle rotor, universal tube holders, and 2-year warranty included in the affordable price. All units ship in 3-5 business days.




### Fixed Angle Rotor

The Apex 4's fixed angle rotor is compatible with recommended IFUs for all common tube manufacturers, including Greiner and BD. Spin four tubes up to 6 mL, or 1.5-2 mL micro BD Microtainers.

Loading Angle

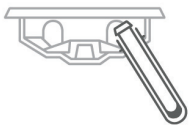

Rotation Angle

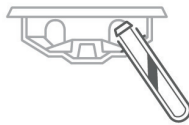

Applications	Max Speed	Max Capacity	Included	Warranty
STAT Chemistry, Coag, PPP	4,400 xg Max RCF 8,500 Max RPM	4 tubes up to 75 mm	Fixed Angle Rotor, Tube Holders	2 years Rotor: lifetime warranty

# DASH APEX 4



## Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	00-053-009-000
Pediatric Tube Adapter, 1.5-2.0mL	03-1-0007-0059 (4 qty)
Stat Rotor Adapter, 13 x 75mm	7713064
1 yr Warranty Extension	03-0-0007-0003
2 yr Warranty Extension	03-0-0007-0004
3 yr Warranty Extension	03-0-0007-0005

Factory Settings			
	RPM	Time	G-Force
1st Setting/Cycle	8,500	30 sec	4,400
2nd Setting/Cycle	8,500	2 min	4,400
3rd Setting/Cycle	8,500	3 min	4,400

Specifications	
Sound level	59 dBA
Dimensions (W x D x H)	7 in x 7 in x 7 in (18 cm x 18 cm x 18 cm)
Weight	4.3 lbs. (1.95 kg)
Power	100 – 240 (+/- 10%) VAC – 50/60 Hz
Max Speed/G-Force	8,500 RPM/4,400 xg



druckerdiagnostics.com  
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)  
+1-814-692-7661