

DASH APEX 6

SET & LOCK STAT CENTRIFUGE

Part # 00-076-009-004



Discover a better way to spin STAT with Drucker's DASH Apex 6 STAT centrifuges. Minimize lab error, reduce turnaround time, and streamline training with the DASH Apex's simple controls and powerful performance.

Easy to Monitor

Prevent forgotten samples and reduce TAT with LED lid lighting. Bright lights indicate when the centrifuge is ready to load (off), running (on), or done (flashing). Lights continue flashing until the lid is opened.

Easy to Use

The simple two-button interface makes it easy to select a preset or program your custom validated STAT cycles. Use all three presets, or set and lock to a single cycle for maximum repeatability. Spin tubes up to 16 x 100 mm or 10 mL with the same rotor and universal tube holders.

Eliminate Wait Time

Further shorten TAT by using a centrifuge array in accessioning to eliminate queuing. The array takes advantage of the DASH Apex 6's compact size to allow samples to be spun as soon as they arrive.

Safe and Durable

Enjoy years of worry-free operation with the DASH Apex 6's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each centrifuge comes ready to run with a horizontal rotor, carriers, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

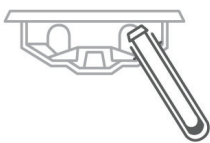

Rotation Angle

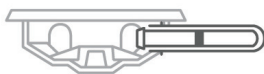

Applications	Max Speed	Max Capacity	Included	Warranty
STAT Chemistry, Coag, PPP	5300 Max RPM 4000 xg Max RCF	6 x up to 16 x 100 mm or 10 mL	Horizontal Rotor 6 x 16 x 100 mm CF Carrier	2 years Extended warranty options available

Designed, built, and supported in the USA

DASH APEX 6



Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles	Time	RPM	G-Force
Cycle 1	3 min	5300	4000
Cycle 2	5 min	5300	4000
Cycle 3	7 min	3800	2000

Part Numbers	
Centrifuge	00-076-009-004
Included	16 x 100 mm Carrier, Carbon Fiber (6 qty) 03-1-0007-0046K
	1-1.5 mL Microtainer Tube Adapter (6 qty) 03-1-0007-0059K
Optional	1 yr Warranty Extension 03-0-0007-0003
	2 yr Warranty Extension 03-0-0007-0004
	3 yr Warranty Extension 03-0-0007-0005

(adapters available for a variety of tube sizes)

Specifications	
Sound Level	61 dBA
Max RPM	5300 RPM
Max G-Force	4000 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30.5 x 35.6 x 22.9 cm
Weight	12 lbs (5.4 kg)
Power	100 - 240 (+/- 10%) - 50/60 Hz



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661