

## DASH FLEX 24 PROGRAMMABLE STAT CENTRIFUGE

Part # 00-184-009-000



Always know when your STAT samples are ready with Drucker's DASH Flex 24 centrifuge. The unique flashing lid lights make it easy to track cycle status and keep TAT down. Customize settings and check every detail on the digital display, and save up to 10 custom cycles to produce high-quality STAT samples.

### Easy to Monitor

LED lid lighting indicates when the centrifuge is ready (off), running (on), or done (flashing), reducing forgotten samples and keeping your average TAT down. The DASH Flex 24's digital display shows cycle parameters and time remaining.

### Easy to Use

Set and save up to 10 programmable time, speed, and g-force settings using the user-friendly interface. The DASH Flex 24 arrives pre-programmed with the most common STAT settings, and can be easily reprogrammed or adjusted on the fly to meet your lab's specimen processing needs.

### High Capacity STAT

The DASH Flex 24 can spin 24 tubes up to 16 x 100 mm or 10 mL, 12 tubes up to 16 x 125 mm or 15 mL, and 6 tubes up to 50 mL using a variety of different carriers and adapters.

### Safe and Durable

Enjoy years of worry-free operation with the DASH's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


### A Great Value

Each DASH centrifuge comes fully equipped and ready to run with a horizontal rotor, carriers, and 2-year warranty included in the affordable price.


### 2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

#### Loading Angle

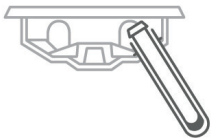

#### Rotation Angle

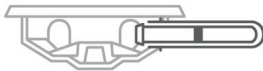

Applications	Max Speed	Max Capacity	Included	Warranty
STAT Chemistry, Coag, PPP	4500 Max RPM 3450 xg Max RCF	24 x up to 16 x 100 mm or 10 mL  12 x up to 16 x 125 mm or 15 mL  6 x 50 mL	Horizontal Rotor  6 x 16 x 125 mm 4-Place CF Carrier	2 years  Extended warranty options available

**Designed, built, and supported in the USA**

# DASH FLEX 24



## Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles	Time	RPM	G-Force
Cycle 1	5 min	4300	3150
Cycle 2	6 min	4300	3150
Cycle 3	7 min	3500	2100

Part Numbers	
Centrifuge	00-184-009-000
Included	16 x 125 mm 4-Place Carrier, Carbon Fiber (6 qty)
	50 mL Carrier, Carbon Fiber (6 qty)
Optional	1-1.5 mL Tube Adapter, Purple (6 qty)
	1 yr Warranty Extension
	2 yr Warranty Extension
	3 yr Warranty Extension

*(adapters available for a variety of tube sizes)*

Specifications	
Sound Level	64 dBA
Max RPM	4500 RPM
Max G-Force	3450 xg
Dimensions (W x D x H)	15 x 17 x 9 in 38.1 x 43.2 x 22.9 cm
Weight	39 lbs (17 kg)
Power	100-240 (+/- 10%) - 50/60 Hz



druckerdiagnostics.com  
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)  
+1-814-692-7661