

## DASH FLEX 6 PROGRAMMABLE STAT CENTRIFUGE

Part # 00-176-009-000

Drucker Diagnostics' DASH Flex 6 introduces the flexible controls and digital display of the best-selling DASH Flex 12 to a smaller body, saving you valuable bench space and allowing you to customize cycle parameters based on your IFUs and unique validation needs.

### Easy to Monitor

Prevent forgotten samples and reduce TAT with LED lid lighting. Bright lights indicate when the centrifuge is ready to load (off), running (on), or done (flashing). Lights continue flashing until the lid is opened. DASH Flex 6 digital display shows cycle parameters and time remaining.

### Easy to Use

Set and save up to 10 programmable time, speed, and g-force settings using the user-friendly interface. The DASH Flex 6 arrives pre-programmed with the most common STAT settings and can spin tubes up to 16 x 100 mm or 10 mL with the same rotor and carriers.

### Eliminate Wait Time

Further shorten TAT by using a centrifuge array in accessioning to eliminate queuing. The array takes advantage of the DASH Flex 6's compact size to allow samples to be spun as soon as they arrive.

### Safe and Durable

Enjoy years of worry-free operation with DASH's brushless motor, reinforced carbon fiber components, and two-year calibration guarantee. Enhance lab safety with imbalance protection and a shatterproof lid, which locks until the centrifuge has stopped.

### A Great Value


Start spinning samples right out of the box, no assembly or additional accessories required. Each centrifuge arrives with a horizontal rotor, carriers, and 2-year warranty included in the affordable price.




### 2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor is compatible with recommended IFUs for all common tube manufacturers, including Greiner and BD.

**Loading Angle**

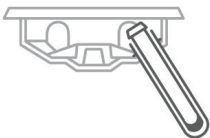

**Rotation Angle**

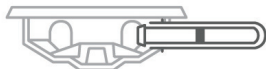

Applications	Max Speed	Max Capacity	Included	Warranty
STAT Chemistry, Coag, PPP	5400 Max RPM 4150 xg Max RCF	6 x up to 16 x 100 mm or 10 mL	Horizontal Rotor  6 x 16 x 100 mm CF Carrier	2 years  Extended warranty options available

# DASH FLEX 6



## Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles	Time	RPM	G-Force
Cycle 1	3 min	5300	4000
Cycle 2	5 min	5300	4000
Cycle 3	10 min	3800	2000

Part Numbers	
Centrifuge	00-176-009-000
Included	16 x 100 mm Carrier, Carbon Fiber (6 qty) 03-1-0007-0046K
	1-1.5 mL Microtainer Tube Adapter (6 qty) 03-1-0007-0059K
Optional	1 yr Warranty Extension 03-0-0007-0003
	2 yr Warranty Extension 03-0-0007-0004
	3 yr Warranty Extension 03-0-0007-0005

(adapters available for a variety of tube sizes)

Specifications	
Sound Level	61 dBA
Max RPM	5400 RPM
Max G-Force	4150 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30.5 x 35.6 x 22.9 cm
Weight	12 lbs (5.4 kg)
Power	100-240 (+/- 10%) - 50/60 Hz



druckerdiagnostics.com  
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)  
+1-814-692-7661