

DASH FLEX 12 PROGRAMMABLE STAT CENTRIFUGE

Part # 00-183-009-000



Always know when your STAT samples are ready with Drucker's DASH Flex 12 centrifuge. The unique flashing lid lights make it easy to track cycle status and keep TAT down. Customize settings and check every detail on the digital display, and save up to 10 custom cycles to produce high-quality STAT samples.

Easy to Monitor

LED lid lighting indicates when the centrifuge is ready (off), running (on), or done (flashing), reducing forgotten samples and keeping your average TAT down. The DASH Flex 12's digital display shows cycle parameters and time remaining.

Easy to Use

Set and save up to 10 programmable time, speed, and g-force settings using the user-friendly interface. The DASH Flex 12 arrives pre-programmed with the most common STAT settings and can spin tubes up to 16 x 100 mm or 10 mL with the same rotor and carriers.

Eliminate Wait Time

The DASH Flex 12 is so compact that it can be used in an array in accessioning. The array will allow spinning samples as soon as they arrive, further shortening your TAT. For a more compact unit, consider the 6-tube DASH Apex 6.

Safe and Durable

Enjoy years of worry-free operation with the DASH's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each DASH centrifuge comes fully equipped and ready to run with a horizontal rotor, carriers, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

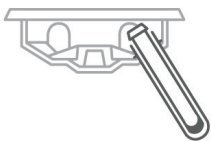

Rotation Angle

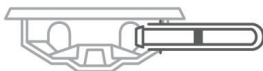

Applications

STAT Chemistry,
Coag, PPP

Max Speed

5200 Max RPM
4000 Max RCF

Max Capacity

12 x up to 16 x 100
mm or 10 mL

Included

Horizontal Rotor
12 x 16 x 100 mm
CF Carrier

Warranty

2 years
Extended warranty
options available

DASH FLEX 12



Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles	Time	RPM	G-Force
Cycle 1	3 min	5200	4000
Cycle 2	5 min	5200	4000
Cycle 3	7 min	3600	2000

Part Numbers	
Centrifuge	00-183-009-000
Included	16 x 100 mm Carrier, Carbon Fiber (12 qty)
	1-1.5 mL Microtainer Tube Adapters (6 qty)
Optional	1 yr Warranty Extension
	2 yr Warranty Extension
	3 yr Warranty Extension

Specifications	
Sound Level	61 dBA
Max RPM	5200 RPM
Max G-Force	4000 xg
Dimensions (W x D x H)	13 x 15 x 9 in 33 x 38.1 x 22.9 cm
Weight	34 lbs (13.7 kg)
Power	100 - 240 (+/- 10%) - 50/60 Hz



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661