

## DASH FLEX 12

### PROGRAMMABLE STAT CENTRIFUGE

Part # 00-183-009-003



Always know when your STAT samples are ready with Drucker's DASH Flex 12 centrifuge. The unique flashing lid lights (patent pending) make it easy to track cycle status and keep TAT down. Customize settings and check every detail on the digital display, and save up to 10 custom cycles to produce high-quality STAT samples.

#### Easy to Monitor

LED lid lighting indicates when the centrifuge is ready (off), running (on), or done (flashing), reducing forgotten samples and keeping your average TAT down. The DASH Flex 12's digital display shows cycle parameters and time remaining.

#### Easy to Use

Set and save up to 10 programmable time, speed, and g-force settings using the user-friendly interface. The DASH Flex 12 arrives pre-programmed with the most common STAT settings and can spin 75-100 mm (3-10 mL) tubes with the same rotor and universal tube holders.

#### Eliminate Wait Time

The DASH Flex 12 is so compact that it can be used in an array in accessioning. The array will allow spinning samples as soon as they arrive, further shortening your TAT. For a more compact unit, consider the 6-tube DASH Apex 6.

#### Safe and Durable

Enjoy years of worry-free operation with the DASH's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


#### A Great Value

Each DASH centrifuge comes fully equipped and ready to run with a horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.


#### 2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

**Loading Angle**

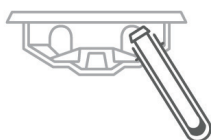

**Rotation Angle**

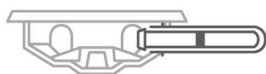

| Applications                 | Max Speed                         | Max Capacity                 | Included                          | Warranty                            |
|------------------------------|-----------------------------------|------------------------------|-----------------------------------|-------------------------------------|
| STAT Chemistry,<br>Coag, PPP | 4,000 xg Max RCF<br>5,200 Max RPM | 12 x<br>75-100 mm<br>3-10 mL | Horizontal Rotor,<br>Tube Holders | 2 years<br>Rotor: lifetime warranty |

**Designed, built, and supported in the USA**

# DASH FLEX 12



## Free 30-Day Evaluation (U.S. Only)

| Part Numbers                               |                        |
|--|------------------------|
| Centrifuge                                 | 00-183-009-003         |
| Replacement 75-100 mm Tube Holders/Buckets | 03-1-0007-0046 (6 qty) |
| Adapters, 0.5 to 1 mL                      | 7713068 (1 qty)        |
| Adapters, 1.5 to 2 mL                      | 7713065 (1 qty)        |
| 1 yr Warranty Extension                    | 03-0-0007-0003         |
| 2 yr Warranty Extension                    | 03-0-0007-0004         |
| 3 yr Warranty Extension                    | 03-0-0007-0005         |

| Factory Settings  |       |       |          |
|-------------------|-------|-------|----------|
|                   | RPM   | Time  | G-Force  |
| 1st Setting/Cycle | 5,200 | 3 min | 4,000 xg |
| 2nd Setting/Cycle | 5,200 | 5 min | 4,000 xg |
| 3rd Setting/Cycle | 3,600 | 7 min | 2,000 xg |

| Specifications         |                                   |
|------------------------|-----------------------------------|
| Sound level            | 61 dBA                            |
| Dimensions (W x D x H) | 13 x 15 x 9 in<br>33 x 38 x 23 cm |
| Weight                 | 30 lbs (13.7 kg)                  |
| Power                  | 95 to 253 VAC – 50/60 Hz          |
| Max Speed/G-Force      | 5200 RPM/4000 xg                  |



druckerdiagnostics.com  
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)  
+1-814-692-7661