

HORIZON 6 FLEX FA FIXED ANGLE ROUTINE CENTRIFUGE

Part # 00-376-009-001

Spin up to 16 x 125 mm or 15 mL tubes and meet fixed angle processing needs for chemistry, urine, and more with HORIZON. Choose the HORIZON 6 Flex FA for its digital display and 10 programmable cycles.



Easy to Use

Program and recall 10 cycles on the HORIZON 6 Flex FA through the digital display. Process tubes up to 16 x 100 mm or 10 mL in the included blue carriers, or tubes up to 16 x 125 mm or 15 mL in the included orange carriers.

Easy to Monitor

Process samples more efficiently with HORIZON's LED lid lighting: off when the centrifuge is ready to load, on when running, and flashing at the end of cycle.

Built to Last

No routine maintenance. Brushless motor, engineered composite components.

Engineered for Safety


Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.

A Great Value


Unbeatable quality for any budget. Each HORIZON centrifuge is ready to run with fixed angle rotor, blue carriers, orange carriers, and 2-year warranty included in the affordable price.

Fixed Angle Rotor

Loading Angle

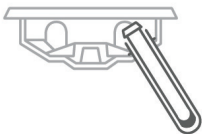

Rotation Angle

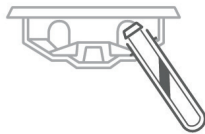

Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag, PPP, Urine	3900 Max RPM 1850 xg Max RCF	6 x up to 16 x 100 mm or 10 mL 6 x up to 16 x 125 mm or 15 mL	Fixed Angle Rotor 6 x 16 x 100 mm Blue Carrier 6 x 15 x 125 mm Orange Carrier	2 years Extended warranty options available

Designed, built, and supported in the USA

HORIZON 6 FLEX FA



*presets based on CLSI guidelines and tube manufacturer IFU

Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles	Time	RPM	With Blue Carrier (Up to 16 x 100 mm or 10 mL tubes)	With Orange Carrier (Up to 16 x 125 mm or 15 mL tubes)
Cycle 1 (Chemistry)	10 min	3900 RPM	1600 xg	1850 xg
Cycle 2 (Coag, PPP)	15 min	3800 RPM	1500 xg	N/A
Cycle 3 (Urine)	5 min	2200 RPM	500 xg	600 xg

Part Numbers	
Centrifuge	00-376-009-001
Included	16 x 100 mm Carrier, Blue (6 qty)
	16 x 125 mm Carrier, Orange (6 qty)
Optional	1-1.5 mL Microtainer Tube Adapter (6 qty)
	1 yr Warranty Extension
	2 yr Warranty Extension
	3 yr Warranty Extension

(adapters available for a variety of tube sizes)

Specifications	
Sound Level	64 dBA
Max RPM	3900 RPM
Max G-Force	1850 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30.5 x 35.6 x 22.9 cm
Weight	12 lbs (5.4 kg)
Power	100 - 240 (+/- 10%) – 50/60 Hz



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661