





Tube Capacity	6 (75, 100, or 125 mm/3-15 mL)
Applications	PRP, Chemistry, PPP, Coag, Urine
Centrifugation	Fixed Angle (max 1,850 xg/3,900 RPM)
Included	Rotor, tube holders
Warranty *	2 years (rotor: lifetime warranty)

* Extended warranty available

HORIZON 6 FLEX FA, HORIZON 6 FA FIXED ANGLE ROUTINE CENTRIFUGES

Spin 125 mm tubes and meet fixed angle processing needs for PRP, urine, and more with HORIZON. Choose the HORIZON 6 Flex FA for its digital display and 10 programmable cycles. Select the HORIZON 6 FA to set and lock your validated settings for maximum standardization.

Easy to Use

Program and recall 10 cycles on the HORIZON 6 Flex FA through the digital display. Set and lock your three preferred cycles on the HORIZON 6 FA through the simple two-button controls. Both spin 75-100 mm (3-10 mL) tubes in the same holder without cushions, or 125 mm (15 mL) tubes in the included additional tube holders.

Easy to Monitor

Process samples more efficiently with HORIZON's LED lid lighting: off when the centrifuge is ready to load, on when running, and flashing at the end of cycle.

Built to Last

No routine maintenance. Brushless motor, engineered composite components.

Engineered for Safety

Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.

A Great Value

Unbeatable quality for any budget. Each HORIZON centrifuge is ready to run with fixed angle rotor, universal 75-100 mm (3-10 mL) tube holders, 125 mm (15 mL) tube holders, and 2-year warranty included in the affordable price.

Choose the HORIZON 6 Flex FA for flexibility & cycle display. Choose the HORIZON 6 FA for simplicity.


Sound Level: 64 dBA

Dimensions (W x D x H): 12 x 14 x 9 in (30 x 36 x 23 cm)

Weight: 12 lbs (5.4 kg)

Power: 95 to 253 V - 50/60 Hz

HORIZON 6 FLEX FA, HORIZON 6 FA

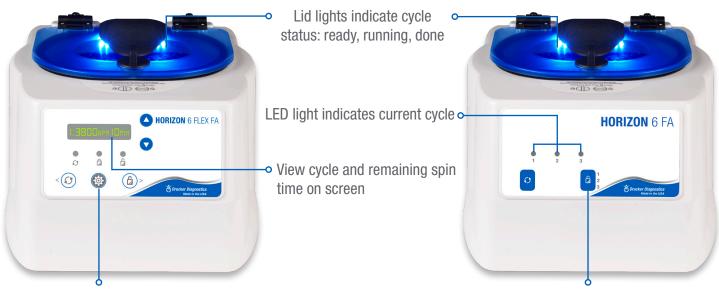

Program your validated cycle, and give your cycle a custom name

Store up to 10 programs, or use default presets for chemistry, coag, and urine*

Select preset cycle: chemistry, coag, urine*
Settings may be locked to a single cycle or
customized at the lab

Free 30-Day Evaluation (U.S. Only)

Part Numbers	HORIZON 6 Flex FA	HORIZON 6 FA
Centrifuge	00-376-009-001	00-276-009-001
Replacement 75-100 mm tube holder/bucket	7713079 (1 qty)	7713079 (1 qty)
Replacement 125 mm orange tube holder/bucket	7713044 (1 qty)	7713044 (1 qty)
SmartView Platform	00-079-009-001	00-079-009-001
1 yr Warranty Extension	03-0-0007-0001	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006	03-0-0007-0006
Custom Settings *	00-100-100-101	00-100-100-101

^{*} Add this part number to your centrifuge order to have your custom validated cycles pre-programmed at the factory

Standard Preset Cycles		With Blue Tube Holder (75-100 mm/3-10 mL tubes)	With Orange Tube Holder (125 mm/15 mL tubes)
1 (Chemistry)	10 min @ 3,900 RPM	1,600 xg	1,850 xg
2 (Coag, PPP)	15 min @ 3,800 RPM	1,500 xg	N/A
3 (PRP, Urine)	5 min @ 2,200 RPM	500 xg	600 xg















druckerdiagnostics.com sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only) +1-814-692-7661



^{*} presets based on CLSI guidelines and tube manufacturer IFU