



## TRUEBOND

### **MULTIPURPOSE VETERINARY CENTRIFUGE**

Part # 00-376-009-005

Using the TrueBond multipurpose veterinary centrifuge, vet laboratories can now spin microhematocrit tubes, fecal, blood, urine, and more in one centrifuge.

#### Spin Any Tube

TrueBond is designed to spin all your tubes, from microhematocrit to urine and everything in between. Includes all necessary adapters.

#### Easy to Use

Minimize training time with the user-friendly interface. The fully digital display makes it easy to customize cycle parameters on the fly, or set and select one of the 10 programmable cycles.

#### Easy to Monitor

Process samples more efficiently with TrueBond's LED lid lighting. Once the unit is started, the display shows remaining spin time.

#### Safe and Durable

Enjoy years of safe, worry-free operation with the brushless motor, imbalance protection, and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


#### A Great Value

Each TrueBond centrifuge is ready to run with a fixed angle rotor, two sets of carriers, microtainer adapters, microhematocrit adapters and reader card, and a 2-year warranty included in the affordable price.



Applications	Max Speed	Max Capacity	Included	Warranty
Blood, Urine, Fecal, Microhematocrit	3900 Max RPM 1850 xg Max RCF	6 x blood, urine, and fecal tubes Up to 16 x 125 mm or 15 mL  8 x Microhematocrits	Fixed Angle Rotor 6 x 16 x 100 mm Blue Carrier 6 x 16 x 125 mm Orange Carrier 6 x 1-1.5 mL Microtainer Adapter 2 x Microhematocrit Adapters Veterinary Microhematocrit Reader Card	2 years Extended warranty options available

# TrueBond

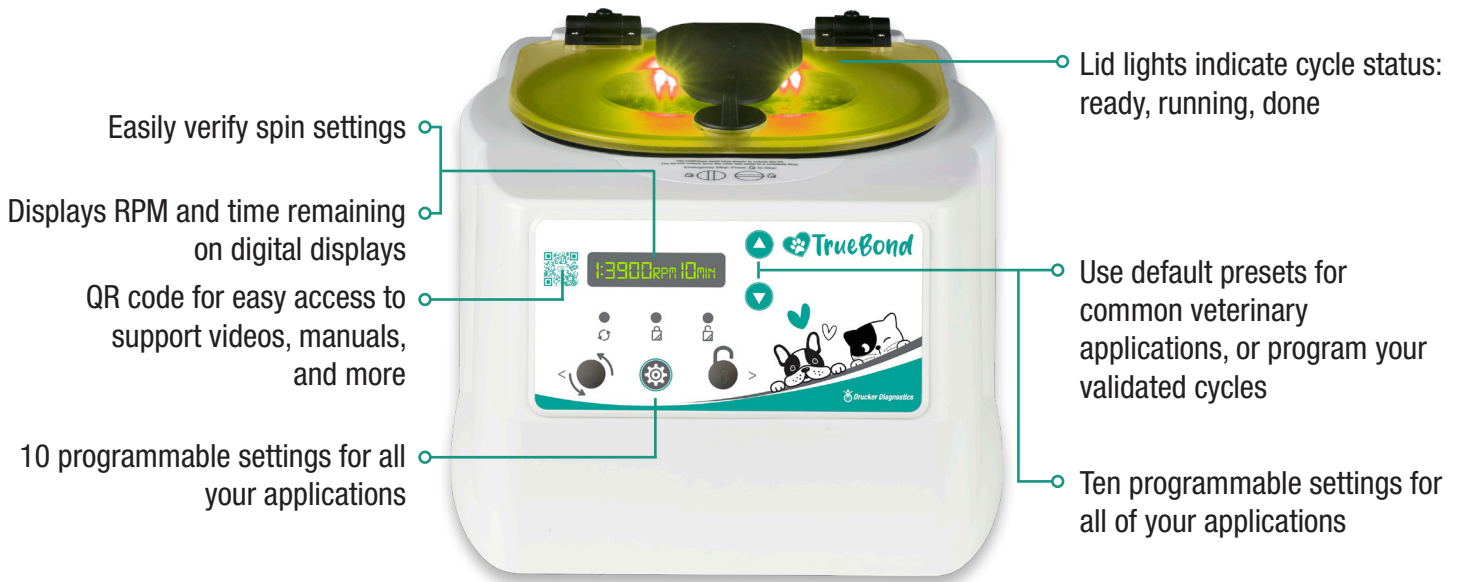

## Free 30-Day Evaluation (U.S. Only)



Standard Preset Cycles	Time	RPM	With Blue Carriers	With Orange Carriers
			(Up to 16 x 100mm or 10 mL tubes)	(Up to 16 x 125 mm or 15 mL tubes)
Cycle 1 (Chemistry)	10 min	3900	1650 xg	N/A
Cycle 2 (Urine Sediment)	5 min	2000	N/A	500 xg
Cycle 3 (Veterinary Fecal)	5 min	1600	300 xg	N/A
Cycle 4 (Microhematocrit)	6 min	3800	1550 xg	N/A

### Part Numbers

Part Numbers		
Centrifuge	00-376-009-005	
Included	16 x 100 mm Carrier, Blue (6 qty)	7713079K
	16 x 125 mm Carrier, Orange (6 qty)	7713044K
	1-1.5 mL Microtainer Tube Adapter (6 qty)	03-1-0007-0059K
	Microhematocrit Adapter (2 qty)	02-004-0-0027K
	Veterinary Microhematocrit Reader Card	02-0006-0-0046
Optional	1 Year Warranty Extension	03-0-0007-0001
	2 Year Warranty Extension	03-0-0007-0002
	3 Year Warranty Extension	03-0-0007-0006

(adapters available for a variety of tube sizes)

### Specifications

Sound Level	64 dBA
Max RPM	3900 RPM
Max G-Force (Blue Carriers)	1600 xg
Max G-Force (Orange Carriers)	1850 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30.5 x 35.6 x 22.9 cm
Weight	12 lbs (5.4 kg)
Power	95 to 253 V – 50/60 Hz



[druckerdiagnostics.com](http://druckerdiagnostics.com)  
[sales@druckerdiagnostics.com](mailto:sales@druckerdiagnostics.com)

+1-866-265-1486 (U.S. only)  
 +1-814-692-7661