



Horizon Centrifuge Model 642B

maximizing bench space and
sample integrity.

OVERVIEW

The Model 642B is our most affordable horizontal centrifuge, but still includes all of the benefits of Drucker's horizontal sample processing. Using our patented horizontal rotor technology, advanced engineering materials, and an advanced cooling system, the 642B ensures that high quality sample separations are achieved & maintained throughout sample transport and analysis.

FEATURES

- **Mechanical Timer:** 1 min. to 30 min allows for a range of applications
- **Horizontal Rotor:** High quality separations for up to six (6) tubes
- **Tube Sizes:** 0.5mL to 10mL; up to 17mm x 100mm
- **Brushless AC Motor:** Years of operation with no routine maintenance
- **Quiet Operation:** A double encased motor and the "stability mount" design provides excellent imbalance stability
- **Cool Operation:** Unique design of the convection based cooling system protects the samples even under the heaviest workloads
- **2 year warrant:** Peace of mind that parts, labor, and support are fully covered


SAFETY

- **Safety Lid Switch:** Centrifuge will not operate with the lid open
- **Positive Lid Latching:** Firm mechanical latching of the lid prevents the escape of debris
- **Clear Shatter Proof Lid:** Allows for safe observation and optical calibrations
- **Regulatory Approvals:** MET marked and certified to UL standard 61010A-1 & 61010A-2-20 and CAN /CSA C22.2 No. 1010.1 & C22.2 No. 1010.2.20
- **Reinforced Side Walls:** Glass fiber reinforced nylon & high-impact polycarbonate provide dual layers of safety protection

APPLICATIONS





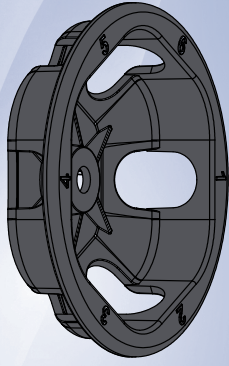

- **Blood Separations:** Chemistry, Coagulation, Serum Separators, and More
- **Urine Separations:** Urinalysis
- **General Purpose Use:** Other Clinical or Analytic Lab Applications

642B CENTRIFUGE ROTOR & ACCESSORIES



6 Place Horizontal Rotor
p/n 7786061

Included with each centrifuge

	QTY Included	Part Number	Comments
	6	7713031	For use with 17 x 100 mm tubes
	6	7713033	For use with 17 x 75 mm tubes
	2	7713066	For use with 13 x 100 mm tubes
	2	7713064	For use with 13 x 75 mm tubes
	2	7713065	For use with tubes between 1.5 mL and 2 mL
	2	7713068	For use with tubes up to 1 mL

Optional Accessories

	1" Cushion Spacer	Available as an accessory	1525	For use with glass vials or tubes
	Tube Holder Cap	Available as an accessory	7713011	For use when additional protection is desired