

HORIZON

HORIZON 6, 12, & 24 *ROUTINE CENTRIFUGES*

Lid lights flash at cycle ends and continue flashing until the lid is opened.

3 cycles preset for Chemistry, Coag, and Urine.

Use presets or program your validated cycles using the two buttons.



Press to select cycle

Press for 5 seconds to lock on a cycle. Repeat to unlock.

Standardize the whole lab on these centrifuges to simplify work flows, training, warranty, and service. Move them around to address bottle necks.

Tubes	Applications	Centrifugation	Included	Warranty
6, 12, or 24 75 & 100 mm	Plasma Serum Urine	Horizontal Max 2,000 xg (RCF)	Rotor Buckets or tube holders	2 years Rotor: lifetime warranty Extended warranty available

Designed, built, and supported in the USA



Horizon 6



Horizon 12




Horizon 24

Part Numbers	Horizon 6	Horizon 12	Horizon 24
Centrifuge	00-276-009-000	00-283-009-000	00-284-009-000
Replacement 75 – 100 mm tube holder/Bucket	7713079	7713079	7713023
50 mL tube holder			7713017
Cap			7713035
SmartView Platform	00-079-009-001		
1 yr Extended Warranty	03-0-0007-0001	03-0-0007-0001	03-0-0007-0001
2 yr Extended Warranty	03-0-0007-0002	03-0-0007-0002	03-0-0007-0002

- Quiet operation
- Positive safety latch
- Maintenance free brushless motor
- Safety lid switch prevents centrifuge from running when lid is open
- Clear shatterproof lid

Cool Operation

Unique design of the motor housing and rotation chamber keeps specimens cool as required by CLSI H21-A5, even under the heaviest workloads.

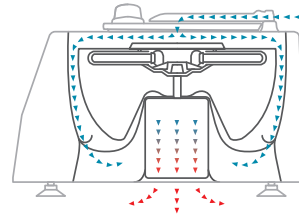

Preset Cycles	Preset Cycles	Horizon 6	Horizon 12	Horizon 24
1 (Chemistry)	10 min @ 1800 xg	3600 RPM	3500 RPM	3300 RPM
2 (Coag – PPP)	15 min @ 1500 xg	3300 RPM	3200 RPM	3000 RPM
3 (Urine)	5 min @ 500 xg	1900 RPM	1800 RPM	1700 RPM

Specifications	Horizon 6	Horizon 12	Horizon 24
Sound level	63.6 dB	55.4 dB	58.2 dB
Dimensions	8 x 11 x 13.5 in 20 x 28 x 34 cm	9 x 12.5 x 14.75 in 23 x 32 x 37 cm	9 x 14.5 x 17 in 23 x 37 x 43 cm
Weight	12 lbs (5.4 kg)	34 lbs. (15 kg)	39 lbs. (17 kg)
Power	100 to 240 V – 50 & 60 Hz		



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-877-231-3115 (U.S. only)
+1-814-692-7661

