

fisherbrand

HORIZON

HORIZON 6, 12, & 24 SET & LOCK ROUTINE CENTRIFUGES

Minimize error, reduce training time, and deliver the best patient outcome with HORIZON. Set and lock your validated cycles to prevent accidental changes. Further reduce potential errors with the option to lock to a single cycle.

Easy to Use

Use the simple, two-button interface to select a preset, or program and record your validated settings on the included label. Process 75-100 mm (3-10 mL) tubes in the same holder without cushions.

Easy to Monitor

Process samples more efficiently with HORIZON's LED lid lighting: off when the centrifuge is ready to load, on when running, and flashing at the end of cycle.

Built to Last

No routine maintenance. Brushless motor, engineered composite components.

Engineered for Safety

Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.

A Great Value

Unbeatable quality for any budget. Each HORIZON centrifuge is ready to run with horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.



Tube Capacity	6, 12, or 24 (75-100 mm/3-10 mL)
Applications	Chemistry, Coag, PPP, PRP, Urine
Centrifugation	Horizontal (max 2,000 xg/3,800 RPM)
Included	Rotor, buckets and tube holders, cycle label
Warranty *	2 years (rotor: lifetime warranty)

* *Extended warranty available*



HORIZON 6




HORIZON 12



HORIZON 24

Designed, built, and supported in the USA

HORIZON 6, 12, & 24

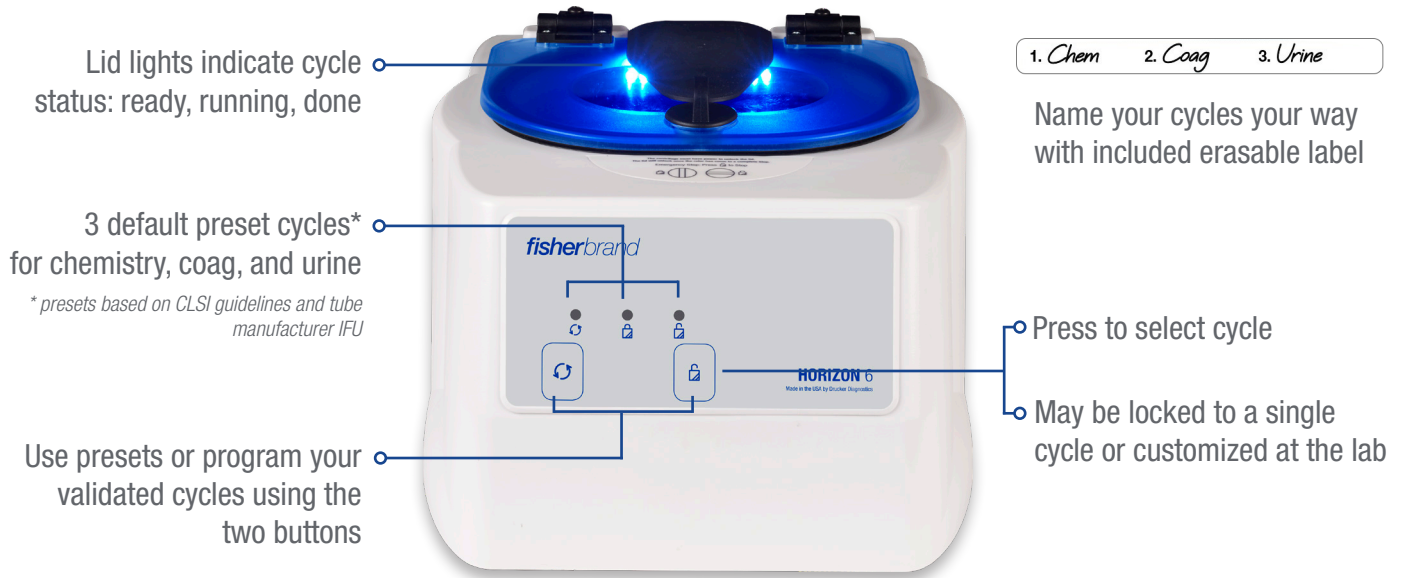

Free 30-Day Evaluation (U.S. Only)

Part Numbers	HORIZON 6	HORIZON 12	HORIZON 24
Centrifuge	FCAT # 22029670	FCAT # 22029674	FCAT # 22029675
Replacement 75-100 mm tube holder/Bucket	FCAT # NC1308605 (1 qty)	FCAT # NC1308605 (1 qty)	FCAT # NC1453932 (1 qty)
50 mL tube holder			FCAT # NC1030468
Cap			7713035
SmartView Platform	FCAT # 22029442		
1 yr Warranty Extension	03-0-0007-0001	03-0-0007-0001	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002	03-0-0007-0002	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006	03-0-0007-0006	03-0-0007-0006
* Custom Settings	00-100-100-101	00-100-100-101	00-100-100-101

* Add this part number to your centrifuge order to have your custom validated cycles pre-programmed at the factory

Standard Preset Cycles		HORIZON 6	HORIZON 12	HORIZON 24
1 (Chemistry)	10 min @ 1,800 xg	3,600 RPM	3,500 RPM	3,300 RPM
2 (Coag, PPP)	15 min @ 1,500 xg	3,300 RPM	3,200 RPM	3,000 RPM
3 (PRP, Urine)	5 min @ 500 xg	1,900 RPM	1,800 RPM	1,700 RPM

Specifications	HORIZON 6	HORIZON 12	HORIZON 24
Sound Level	64 dBA	56 dBA	59 dBA
Dimensions (W x D x H)	12 x 14 x 9 in 30 x 36 x 23 cm	13 x 15 x 9 in 33 x 38 x 23 cm	15 x 17 x 9 in 38 x 43 x 23 cm
Weight	12 lbs (5.4 kg)	30 lbs (13.7 kg)	39 lbs (17 kg)
Power	95 to 253 V – 50/60 Hz		

FDA LISTED



RoHS
Compliant



ISO13485
certified



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661

