

614B

FIXED ANGLE SINGLE-SPEED CENTRIFUGE

Part # 00-004-009-001

With a single g-force setting and adjustable mechanical timer, the 614B provides runtime flexibility and fixed angle processing at an affordable price. Use anywhere specimens are analyzed on site.

A Great Value

Each centrifuge arrives ready to run with fixed angle rotor, tube holders, and 2-year warranty included in the affordable price.

Easy to Use

Set right every time, no need to select from a variety of complicated options. The 614B offers a single g-force setting of 1,100 xg and an easily adjustable mechanical dial with runtimes from 1 to 30 minutes.

Cool Operation

Unique design of the motor housing and rotation chamber keeps specimens cool during processing.


Built to Last

Enjoy years of maintenance-free operation with the 614B's reinforced components and brushless motor.


Fixed-Angle Rotor Included

Used for spinning 75 mm, 100 mm, and 125 mm tubes in specially designed tube holders.

Loading Angle

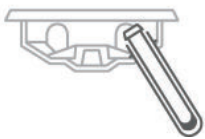

Rotation Angle

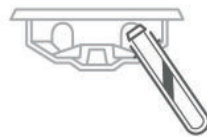

Fixed Angle vs Horizontal Centrifugation

Fixed-angle centrifuges produce gel layers that may be susceptible to rupture during transport. These damaged samples must be re-centrifuged – a process the CLSI does not recommend due to its known negative effects on analytical results. For better gel layer formation, especially for tubes that will be transported, please consider horizontal centrifugation with the **Model 642E** or **HORIZON** series.


Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag	1,200 xg Max RCF 3,250 Max RPM	6 x 75-125 mm 3-15 mL	Fixed Angle Rotor, Tube Holders	2 years Rotor: lifetime warranty

Designed, built, and supported in the USA

614B Centrifuge

Part Numbers	614B
Centrifuge	00-004-009-001
75-100 mm tube holders	7713079K (6 qty)
125 mm tube holders	7713032K (6 qty)
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Specifications	614B
Sound level	55 dBA
Dimensions (W x D x H)	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	11 lbs (5 kg)
Max RPM	3250 RPM
Max G-Force	1200 xg
Power	115 V – 50/60 Hz

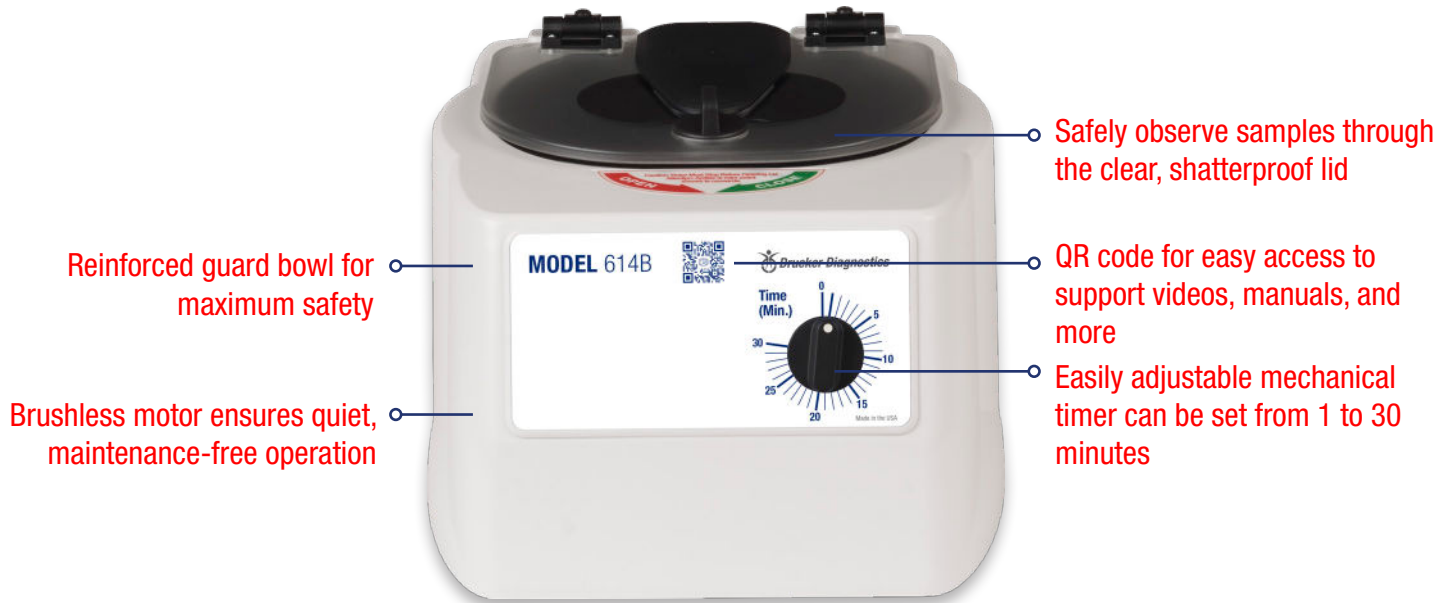

OTHER DRUCKER MODELS



642E

World's #1 clinical lab centrifuge
For blood



HORIZON 6 FA

3 preset cycles, 6 tube capacity
For blood and urine



HORIZON 6 Flex FA

Digital display, 10 cycles, 6 tube capacity
For blood and urine



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661

