

DASH FLEX 24

PROGRAMMABLE STAT CENTRIFUGE

Part # 00-184-009-000



Always know when your STAT samples are ready with Drucker's DASH Flex 24 centrifuge. The unique flashing lid lights (patent pending) make it easy to track cycle status and keep TAT down. Customize settings and check every detail on the digital display, and save up to 10 custom cycles to produce high-quality STAT samples.

Easy to Monitor

LED lid lighting indicates when the centrifuge is ready (off), running (on), or done (flashing), reducing forgotten samples and keeping your average TAT down. The DASH Flex 24's digital display shows cycle parameters and time remaining.

Easy to Use

Set and save up to 10 programmable time, speed, and g-force settings using the user-friendly interface. The DASH Flex 24 arrives pre-programmed with the most common STAT settings, and can be easily reprogrammed or adjusted on the fly to meet your lab's specimen processing needs.

High Capacity STAT

The DASH Flex 24 can spin standard 75-100 mm (3-10 mL) blood tubes, 125 mm (15 mL) long tubes, and tubes with volumes up to 50 mL using a variety of different buckets - no small adapters!

Safe and Durable

Enjoy years of worry-free operation with the DASH's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each DASH centrifuge comes fully equipped and ready to run with a horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

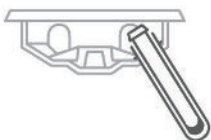

Rotation Angle

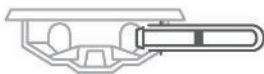

Applications

STAT Chemistry,
Coag, PPP

Max Speed

3,450 xg Max RCF
4,500 Max RPM

Max Capacity

75-100 mm/3-10 mL
24 x Standard Tubes
12 x 125 mm/15 mL
6 x 50 mL


Included

Horizontal Rotor,
Buckets, Tube
Holders

Warranty

2 years
Rotor: lifetime warranty

DASH FLEX 24

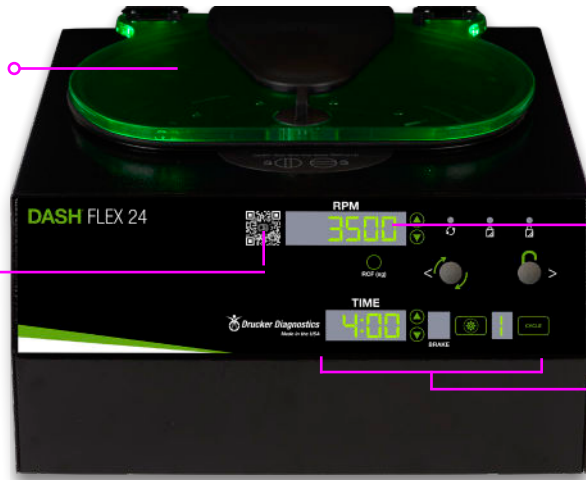

LED lid lights indicate cycle status: ready, running, done

QR code for easy access to support videos, manuals, and more

View RPM and remaining spin time on the digital display

Program and name your custom validated cycles

Store up to 10 settings

Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	00-184-009-000
Replacement Tube Holders/Buckets	02-004-0-0012K (pack of 6)
50 mL Tube Holders/Buckets	02-004-0-0014K (pack of 6)
Adapters, 1.5 to 2 mL	7713065 (1 qty)
1 yr Warranty Extension	03-0-0007-0003
2 yr Warranty Extension	03-0-0007-0004
3 yr Warranty Extension	03-0-0007-0005

Factory Settings			
	RPM	Time	G-Force
1st Setting/Cycle	4,300	5 min	3,000 xg
2nd Setting/Cycle	4,300	6 min	3,000 xg
3rd Setting/Cycle	3,500	7 min	2,000 xg

Specifications	
Sound level	61 dBA
Dimensions (W x D x H)	15 x 17 x 9 in 38 x 43 x 23 cm
Weight	31 lbs (14 kg)
Power	95 to 253 VAC – 50/60 Hz
Max Speed/G-Force	4500 RPM/3450 xg



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661

