

## HORIZON 6 FA

### FIXED ANGLE ROUTINE CENTRIFUGE

Part # 00-276-009-001

Spin 125 mm tubes and meet fixed angle processing needs for onsite chemistry, urine, and more with HORIZON. Choose the HORIZON 6 FA to set and lock your validated settings for maximum standardization.

#### Easy to Use

Use the simple, two-button interface to select a preset, or program and record your validated settings on the included label. Spin 75-100 mm (3-10 mL) tubes in the same holder without cushions, or 125 mm (15 mL) tubes in the included additional tube holders.

#### Easy to Monitor

Process samples more efficiently with HORIZON 6 FA's LED lid lighting: off when the centrifuge is ready to load, on when running, and flashing at the end of cycle.

#### Built to Last

No routine maintenance. Brushless motor, engineered composite components.

#### Engineered for Safety

Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.


#### A Great Value

Unbeatable quality for any budget. Each HORIZON 6 FA centrifuge is ready to run with fixed angle rotor, universal 75-100 mm (3-10 mL) tube holders, 125 mm (15 mL) tube holders, and 2-year warranty included in the affordable price.




#### Fixed Angle Rotor

##### Loading Angle

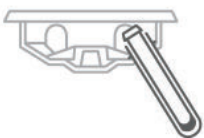


##### Rotation Angle

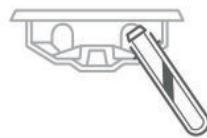

Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag, PPP, Urine	1,850 xg Max RCF 3,900 Max RPM	6 x 75-125 mm 3-15 mL	Fixed Angle Rotor, Tube Holders, Cycle Label	2 years Rotor: lifetime warranty

**Designed, built, and supported in the USA**

# HORIZON 6 FA

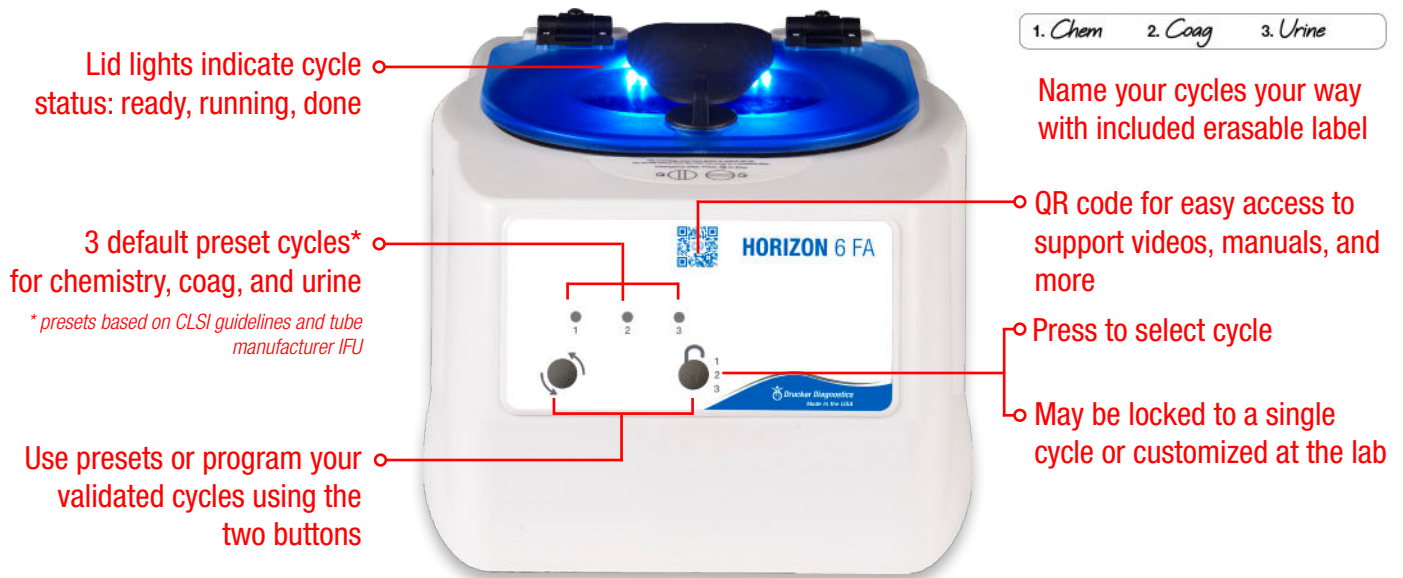

## Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	00-276-009-001
Replacement 75-100 mm Tube Holders	7713079K (pack of 6)
Replacement 125 mm Tube Holders	7713044K (pack of 6)
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Standard Preset Cycles		With Blue Tube Holder (75-100 mm/3-10 mL tubes)	With Orange Tube Holder (125 mm/15 mL tubes)
1 (Chemistry)	10 min @ 3,900 RPM	1,600 xg	1,850 xg
2 (Coag, PPP)	15 min @ 3,800 RPM	1,500 xg	N/A
3 (Urine)	5 min @ 2,200 RPM	500 xg	600 xg

Specifications	
Sound Level	64 dBA
Max RPM/G-Force	3900 RPM/1600 xg/1850 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	12 lbs (5.4 kg)
Power	95 to 253 V – 50/60 Hz



druckerdiagnostics.com  
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)  
+1-814-692-7661

