

HORIZON 12

SET & LOCK ROUTINE CENTRIFUGE

Part # 00-283-009-000



Minimize error, reduce training time, and deliver the best patient outcome with HORIZON 12. Set and lock your validated cycles to prevent accidental changes. Further reduce potential errors with the option to lock to a single cycle.

Easy to Use

Use the simple, two-button interface to select a preset, or program and record your validated settings on the included label. Process 75-100 mm (3-10 mL) tubes in the same holder without cushions.

Easy to Monitor

Process samples more efficiently with HORIZON's LED lid lighting: off when the centrifuge is ready to load, on when running, and flashing at the end of cycle.

Built to Last

No routine maintenance. Brushless motor, engineered composite components.

Engineered for Safety

Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.


A Great Value

Unbeatable quality for any budget. Each HORIZON 12 centrifuge is ready to run with horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

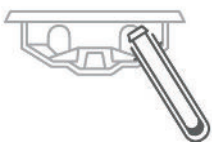

Rotation Angle

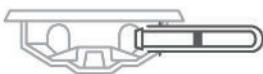

| Applications | Max Speed | Max Capacity | Included | Warranty |
|-----------------------------|-----------------------------------|------------------------------|-----------------------------------|-------------------------------------|
| Chemistry, Coag, PPP, Urine | 2,000 xg Max RCF 3,700 Max RPM | 12 x 75-100 mm 3-10 mL | Horizontal Rotor, Tube Holders | 2 years Rotor: lifetime warranty |

Designed, built, and supported in the USA

HORIZON 12



Free 30-Day Evaluation (U.S. Only)

| Part Numbers | |
|--|----------------------|
| Centrifuge | 00-283-009-000 |
| Replacement 75-100 mm Tube Holders/Buckets | 7713079K (pack of 6) |
| 1 yr Warranty Extension | 03-0-0007-0001 |
| 2 yr Warranty Extension | 03-0-0007-0002 |
| 3 yr Warranty Extension | 03-0-0007-0006 |

| Standard Preset Cycles | Time | G-Force | RPM |
|------------------------|--------|----------|-----------|
| 1 (Chemistry) | 10 min | 1,800 xg | 3,500 RPM |
| 2 (Coag, PPP) | 15 min | 1,500 xg | 3,200 RPM |
| 3 (Urine) | 5 min | 500 xg | 1,800 RPM |

| Specifications | |
|------------------------|-----------------------------------|
| Sound Level | 56 dBA |
| Max RPM/G-Force | 3700 RPM/2000 xg |
| Dimensions (W x D x H) | 13 x 15 x 9 in 33 x 38 x 23 cm |
| Weight | 30 lbs (13.7 kg) |
| Power | 95 to 253 V – 50/60 Hz |



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661

