



Centrifuge Product Catalog



VWR One Spin Clinical Centrifuge

VWR One Spin Mobile Clinical Centrifuge

VWR One Spin Clinical Centrifuge with Timer

VWR Flexible 6 Clinical Centrifuge

VWR Flexible 6 Fixed Angle Clinical Centrifuge

VWR Flexible 12 Clinical Centrifuge

VWR Flexible 24 Clinical Centrifuge


VWR STAT 6 Clinical Centrifuge

VWR STAT 12 Clinical Centrifuge

VWR STAT Coag Clinical Centrifuge

VWR SEROLOGY 12 Clinical Centrifuge

Guide to Selecting a Centrifuge

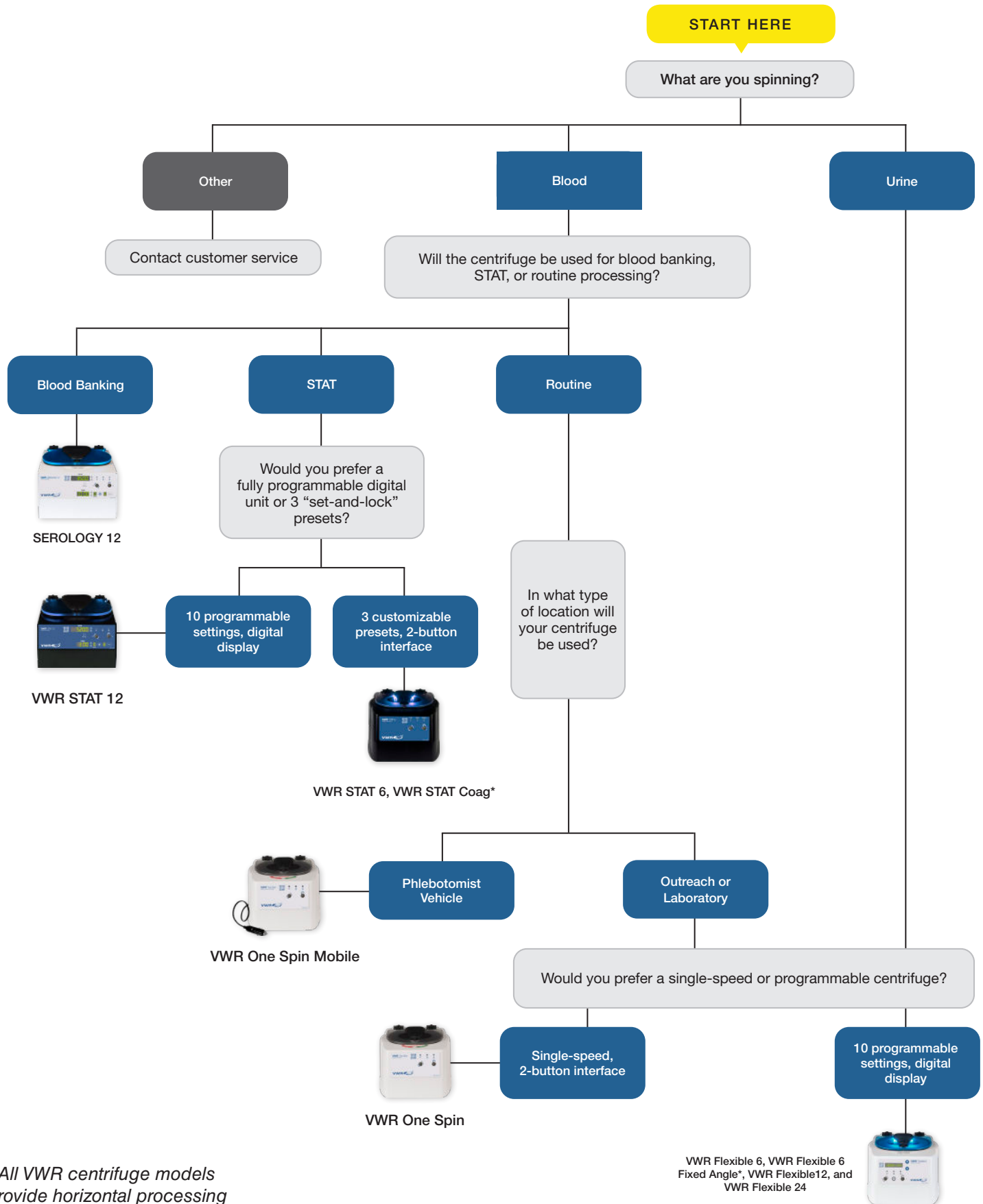

* All VWR centrifuge models provide horizontal processing unless otherwise specified

VWR Flexible 6, VWR Flexible 6 Fixed Angle*, VWR Flexible 12, and VWR Flexible 24

Request a Free 30 Day Evaluation (U.S. Only)

CENTRIFUGE MODEL	PART NUMBER
VWR STAT Coag Clinical Centrifuge*	76795-588
VWR STAT 12 Clinical Centrifuge	76795-576
VWR STAT 6 Clinical Centrifuge	76795-572
VWR Flexible 6 Clinical Centrifuge	76795-578
VWR Flexible 6 Fixed Angle Clinical Centrifuge*	76795-580
VWR Flexible 12 Clinical Centrifuge	76795-582
VWR Flexible 24 Clinical Centrifuge	76795-584
VWR SEROLOGY 12 Clinical Centrifuge	76795-586
VWR One Spin Clinical Centrifuge	76795-570
VWR One Spin Mobile Clinical Centrifuge	76795-574
VWR One Spin Clinical Centrifuge with Timer	76795-568

**fixed angle processing*



vwr.com | sales@vwr.com | 1-800-932-5000



Protected by U.S. Patents #6.811.531. #D718,463S #D734,489S Other Patents Pending.



ONE SPIN

VWR® One Spin Clinical Centrifuge *SINGLE-SPEED ELECTRONIC CENTRIFUGE*

Part # 76795-570

The VWR® One Spin Clinical Centrifuge is the **world's best-selling clinical laboratory centrifuge**. Easy two-button operation ensures repeatability and standardization, while advanced design delivers high-quality blood separation every time. Its preset speed and runtime make the VWR® One Spin Clinical Centrifuge the ideal choice for consistency in doctor's offices, patient service centers, and outreach locations.

Easy to Use

Spin samples with the press of just one button. The VWR® One Spin Clinical Centrifuge processes 75-100 mm (3-10 mL) tubes for chemistry, coag, and PPP.

Simple Settings

The VWR® One Spin Clinical Centrifuge features a preset runtime of 10 minutes at 1,600 xg for all common chemistry needs.

Built to Last

No routine maintenance. The brushless motor and unique convection-based cooling system ensure years of high-quality separations.

Engineered for Safety

Protect your lab with the clear, shatterproof lid, which locks until the centrifuge has stopped completely.

A Great Value


Get the world's #1 clinical laboratory centrifuge at a price that fits any budget. The unit arrives ready for use with the horizontal rotor, tube holders, and 2-year warranty included in the affordable price.




2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

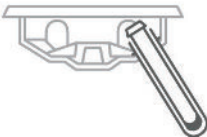

Rotation Angle

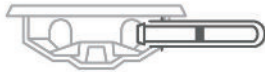

Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag, PPP	1,600 xg Max RCF 3,380 Max RPM	6 x 75-100 mm 3-10 mL	Horizontal Rotor, Tube Holders	2 years Rotor: lifetime warranty

Designed, built, and supported in the USA

VWR® One Spin Clinical Centrifuge

Part Numbers	
Centrifuge	76795-570
75-100 mm tube holders	7713079K (6 qty)
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Specifications	
Sound level	64 dBA
Dimensions	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	11 lbs (5 kg)
Max RPM	3380 RPM
Max G-Force	1600 xg
Power	115 V – 50/60 Hz



Safety lid automatically locks during operation and will not open until the centrifuge has stopped

QR code for easy access to support videos, manuals, and more

Brushless motor ensures quiet, maintenance-free operation

Lights illuminate to indicate running, locked, unlocked

Simple two-button controls

Reinforced guard bowl for maximum safety

OTHER MODELS



VWR® One Spin Mobile Clinical Centrifuge

For mobile centrifugation
For use with 12V outlet



VWR® One Spin Clinical Centrifuge With Timer

Adjustable time, fixed speed
For spinning blood



VWR® Flexible 6 Clinical Centrifuge

Digital display, available for 6, 12, 24 tubes
10 cycles, for blood and urine



vwr.com | sales@vwr.com | 1-800-932-5000



ONE SPIN

VWR® One Spin Mobile Clinical Centrifuge 12V DC MOBILE CENTRIFUGE

Part # 76795-574



The VWR® One Spin Mobile Clinical Centrifuge facilitates efficient routes and longer range for mobile phlebotomists by providing centrifugation in the vehicle. Spin specimens within two hours of draw and run the most convenient route without the need to return to the central lab. Powered directly from a vehicle's 12V electrical supply, the VWR® One Spin Mobile Clinical Centrifuge provides horizontal centrifugation for high quality gel tube separation and two-button operation for maximum repeatability.

Easy to Use

The simple interface runs the preset cycle (10 minutes at 3,400 RPM/1,600 xg) at the touch of a button. Runtime can be adjusted from 1 to 30 minutes.

Designed for Mobile Centrifugation

With in-vehicle centrifugation, phlebotomists are free to plan a day-long route without returning to the lab. This affordable device can greatly increase range of service and productivity without sacrificing specimen quality or patient outcome.

Built to Last

No routine maintenance. The brushless motor and unique convection-based cooling system ensure years of high-quality separations.

Engineered for Safety

Ensure maximum safety with the clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each unit arrives ready to run with horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out rotor design incorporates a unique test tube holder that produces horizontally separated samples, ideal for transportation.

Loading Angle

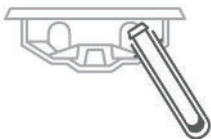

Rotation Angle

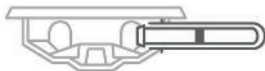

Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag, PPP	1,650 xg Max RCF 3,400 Max RPM	6 x 75-100 mm 3-10 mL	Horizontal Rotor, Tube Holders	2 years Rotor: lifetime warranty

Designed, built, and supported in the USA

VWR® One Spin Mobile Clinical Centrifuge

Part Numbers	
Centrifuge	76795-574
75-100 mm tube holders	7713079K (6 qty)
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Specifications	
Sound level	65 dBA
Dimensions	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	12 lbs (5.4 kg)
Max RPM	3400 RPM
Max G-Force	1650 xg
Power	12V



OTHER MODELS



VWR® One Spin Clinical Centrifuge With Timer

Adjustable time, fixed speed
For spinning blood



VWR® One Spin Clinical Centrifuge

Fixed time, fixed speed
For spinning blood



VWR® Flexible 6 Clinical Centrifuge

Digital display, available for 6, 12, 24 tubes
10 cycles, for blood and urine



vwr.com | sales@vwr.com | 1-800-932-5000

VWR® One Spin Clinical Centrifuge with Timer *HORIZONTAL SINGLE-SPEED CENTRIFUGE*

Part # 76795-568

Produce high-quality sample separations at our most affordable price. The VWR® One Spin Clinical Centrifuge with Timer's simple, single-speed settings and adjustable mechanical timer meet a variety of blood specimen preparation needs. The unique rotor design allows for horizontal processing with minimal maintenance and no hassle.

A Great Value

Unbeatable quality for any budget. Each unit arrives ready to run with a rotor, tube holders, and 2-year warranty included in the affordable price.

Easy to Use

The VWR® One Spin Clinical Centrifuge with Timer offers powerful single-speed processing with an easily adjustable mechanical timer. Set to run at a constant 1,600 xg, the 642B allows selection of cycle time from 1 to 30 minutes.

Cool Operation

Unique design of the motor housing and rotation chamber keeps specimens cool during processing.

Built to Last

Reinforced components and brushless motor provide years of maintenance-free operation. The clear, shatterproof lid enhances lab safety.

For Better Repeatability


When consistency is key, please consider the single-speed, electronic **VWR® One Spin Clinical Centrifuge**. The VWR® One Spin Clinical Centrifuge is the world's best-selling clinical laboratory centrifuge and features a simple, two-button interface and preset runtime for all common chemistry needs.




2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

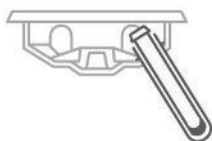

Rotation Angle

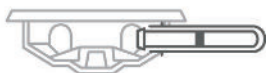

Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag, PPP	1,600 xg Max RCF 3,380 Max RPM	6 x 75-100 mm 3-10 mL	Horizontal Rotor, Tube Holders	2 years Rotor: lifetime warranty

VWR® One Spin Clinical Centrifuge with Timer

Part Numbers	
Centrifuge	76795-568
75-100 mm tube holders	7713079K (6 qty)
125mm tube holders	7713032K (6 qty)
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Specifications	
Sound level	55 dBA
Dimensions	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	11 lbs (5 kg)
Max RPM	3250 RPM
Max G-Force	1200 xg
Power	115 V – 50/60 Hz

Safely observe samples through the clear, shatterproof lid

Reinforced guard bowl and brushless motor for quiet, maintenance-free operation



QR code for easy access to support videos, manuals, and more

Easily adjustable mechanical timer can be set from 1 to 30 minutes

Consistent processing at 1,600 xg (3,400 RPM)

OTHER MODELS



VWR® One Spin Mobile Clinical Centrifuge

For mobile centrifugation
For use with 12V outlet



VWR® One Spin Clinical Centrifuge

Fixed time, fixed speed
For spinning blood



VWR® Flexible 6 Clinical Centrifuge

Digital display, available for 6, 12, 24 tubes
10 cycles, for blood and urine



vwr.com | sales@vwr.com | 1-800-932-5000

VWR® FLEXIBLE 6 CLINICAL CENTRIFUGE PROGRAMMABLE ROUTINE CENTRIFUGE

Part # 76795-578

Satisfy the most diverse processing requirements with the VWR® Flexible 6 Clinical Centrifuge. Customize settings and check every detail on the digital display, or pre-program up to 10 cycles.

Easy to Use

The fully digital display makes it easy to customize cycle parameters on the fly, or set and select one of the 10 programmable cycles. Spin 75-100 mm (3-10 mL) tubes without switching tube holders or adding cushions.

Easy to Monitor

Process samples more efficiently with VWR® Flexible's LED lid lighting: off when the centrifuge is ready to load, on when running, and flashing at the end of cycle. Once the unit is started, the display shows remaining spin time.

Built to Last

No routine maintenance. Brushless motor, engineered composite components.

Engineered for Safety

Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.

A Great Value


Unbeatable quality for any budget. Each VWR® Flexible 6 Clinical Centrifuge is ready to run with horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.




2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

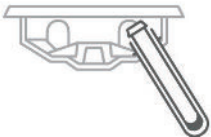

Rotation Angle

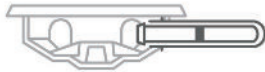

Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag, PPP, Urine	2,000 xg Max RCF 3,800 Max RPM	6 x 75-100 mm 3-10 mL	Horizontal Rotor, Tube Holders	2 years Rotor: lifetime warranty

Designed, built, and supported in the USA

VWR® FLEXIBLE 6 CLINICAL CENTRIFUGE



Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	76795-578
Replacement 75-100 mm Tube Holders/Buckets	7713079K (pack of 6)
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Standard Preset Cycles	Time	G-Force	RPM
1 (Chemistry)	10 min	1,800 xg	3,600 RPM
2 (Coag, PPP)	15 min	1,500 xg	3,300 RPM
3 (Urine)	5 min	500 xg	1,900 RPM

Specifications	
Sound Level	64 dBA
Max RPM/G-Force	3800 RPM/2000 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	12 lbs (5.4 kg)
Power	95 to 253 V – 50/60 Hz



vwr.com | sales@vwr.com | 1-800-932-5000

VWR® Flexible 6 Fixed Angle Clinical Centrifuge *FIXED ANGLE ROUTINE CENTRIFUGE*

Part # 76795-580

Spin 125 mm tubes and meet fixed angle processing needs for chemistry, urine, and more with the VWR® Flexible 6 Fixed Angle Clinical Centrifuge. Choose the VWR® Flexible 6 Fixed Angle Clinical Centrifuge for its digital display and 10 programmable cycles.



Easy to Use

Program and recall 10 cycles on the VWR® Flexible 6 Fixed Angle Clinical Centrifuge through the digital display. Spin 75-100 mm (3-10 mL) tubes in the same holder without cushions, or 125 mm (15 mL) tubes in the included additional tube holders.

Easy to Monitor

Process samples more efficiently with VWR® Flexible's LED lid lighting: off when the centrifuge is ready to load, on when running, and flashing at the end of cycle.

Built to Last

No routine maintenance. Brushless motor, engineered composite components.

Engineered for Safety


Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.

A Great Value


Unbeatable quality for any budget. Each VWR® Flexible 6 Fixed Angle Clinical Centrifuge is ready to run with fixed angle rotor, universal 75-100 mm (3-10 mL) tube holders, 125 mm (15 mL) tube holders, and 2-year warranty included in the affordable price.

Fixed Angle Rotor

Loading Angle

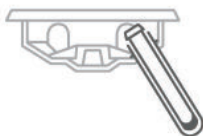

Rotation Angle

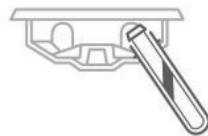

Applications

Chemistry, Coag, PPP, Urine

Max Speed

1,850 xg Max RCF
3,900 Max RPM

Max Capacity

6 x
75-125 mm
3-15 mL

Included

Fixed Angle Rotor,
Tube Holders

Warranty

2 years
Rotor: lifetime warranty

VWR® Flexible 6 Fixed Angle Clinical Centrifuge



Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	76795-580
Replacement 75-100 mm Tube Holders	7713079K (pack of 6)
Replacement 125 mm Tube Holders	7713044K (pack of 6)
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Standard Preset Cycles		With Blue Tube Holder (75-100 mm/3-10 mL tubes)	With Orange Tube Holder (125 mm/15 mL tubes)
1 (Chemistry)	10 min @ 3,900 RPM	1,600 xg	1,850 xg
2 (Coag, PPP)	15 min @ 3,800 RPM	1,500 xg	N/A
3 (Urine)	5 min @ 2,200 RPM	500 xg	600 xg

Specifications	
Sound Level	64 dBA
Max RPM/G-Force	3900 RPM/1600 xg/1850 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	12 lbs (5.4 kg)
Power	95 to 253 V – 50/60 Hz



vwr.com | sales@vwr.com | 1-800-932-5000



FLEXIBLE

VWR® FLEXIBLE 12 CLINICAL CENTRIFUGE PROGRAMMABLE ROUTINE CENTRIFUGE

Part # 76795-582

Satisfy the most diverse processing requirements with the VWR® Flexible 12 Clinical Centrifuge. Customize settings and check every detail on the digital display, or pre-program up to 10 cycles.



Easy to Use

The fully digital display makes it easy to customize cycle parameters on the fly, or set and select one of the 10 programmable cycles. Spin 75-100 mm (3-10 mL) tubes without switching tube holders or adding cushions.

Easy to Monitor

Process samples more efficiently with VWR® Flexible's LED lid lighting: off when the centrifuge is ready to load, on when running, and flashing at the end of cycle. Once the unit is started, the display shows remaining spin time.

Built to Last

No routine maintenance. Brushless motor, engineered composite components.

Engineered for Safety

Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.


A Great Value

Unbeatable quality for any budget. Each VWR® Flexible 12 Clinical Centrifuge is ready to run with horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

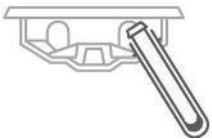

Rotation Angle

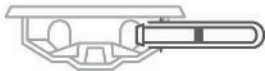

Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag, PPP, Urine	2,000 xg Max RCF 3,700 Max RPM	12 x 75-100 mm 3-10 mL	Horizontal Rotor, Tube Holders	2 years Rotor: lifetime warranty

Designed, built, and supported in the USA

VWR® FLEXIBLE 12 CLINICAL CENTRIFUGE



Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	76795-582
Replacement 75-100 mm Tube Holders/Buckets	7713079K (pack of 6)
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Standard Preset Cycles	Time	G-Force	RPM
1 (Chemistry)	10 min	1,800 xg	3,500 RPM
2 (Coag, PPP)	15 min	1,500 xg	3,200 RPM
3 (Urine)	5 min	500 xg	1,800 RPM

Specifications	
Sound Level	56 dBA
Max RPM/G-Force	3700 RPM/2000 xg
Dimensions (W x D x H)	13 x 15 x 9 in 33 x 38 x 23 cm
Weight	30 lbs (13.7 kg)
Power	95 to 253 V – 50/60 Hz



vwr.com | sales@vwr.com | 1-800-932-5000



FLEXIBLE

VWR® FLEXIBLE 24 CLINICAL CENTRIFUGE PROGRAMMABLE ROUTINE CENTRIFUGE

Part # 76795-584

Satisfy the most diverse processing requirements with the VWR® Flexible 24 Clinical Centrifuge. Customize settings and check every detail on the digital display, or pre-program up to 10 cycles.



Easy to Use

The fully digital display makes it easy to customize cycle parameters on the fly, or set and select one of the 10 programmable cycles. Process 75-100 mm (3-10 mL) tubes and 125 mm (15 mL) in the same holder without cushions. Simply swap holders to run 50 mL tubes.

Easy to Monitor

Process samples more efficiently with VWR® Flexible's LED lid lighting: off when the centrifuge is ready to load, on when running, and flashing at the end of cycle. Once the unit is started, the display shows remaining spin time.

Built to Last

No routine maintenance. Brushless motor, engineered composite components.

Engineered for Safety

Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.


A Great Value

Unbeatable quality for any budget. Each VWR® Flexible 24 Clinical Centrifuge is ready to run with horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

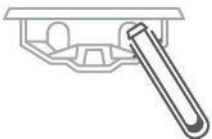

Rotation Angle

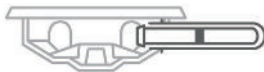

Applications	Max Speed	Max Capacity	Included	Warranty
Chemistry, Coag, PPP, Urine	2,000 xg Max RCF 3,400 Max RPM	75-100 mm/3-10 mL 24 x Standard Tubes 12 x 125 mm/15 mL 6 x 50 mL	Horizontal Rotor, Tube Holders	2 years Rotor: lifetime warranty

Designed, built, and supported in the USA

VWR® FLEXIBLE 24 CLINICAL CENTRIFUGE



LED lid lights indicate cycle status: ready, running, done

QR code for easy access to support videos, manuals, and more

View RPM and remaining spin time on the digital display

Program and name your custom validated cycles

Store up to 10 settings

Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	76795-584
Replacement 75-100 mm Tube Holders/Buckets	7713079K (pack of 6)
50 mL tube holders	7713037K (pack of 6)
Cap	7713035K (pack of 6)
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Standard Preset Cycles	Time	G-Force	RPM
1 (Chemistry)	10 min	1,800 xg	3,300 RPM
2 (Coag, PPP)	15 min	1,500 xg	3,000 RPM
3 (Urine)	5 min	500 xg	1,700 RPM

Specifications	
Sound Level	59 dBA
Max RPM/G-Force	3400 RPM/2000 xg
Dimensions (W x D x H)	15 x 17 x 9 in 38 x 43 x 23 cm
Weight	31 lbs (14 kg)
Power	95 to 253 V – 50/60 Hz



vwr.com | sales@vwr.com | 1-800-932-5000

VWR® STAT 6 Clinical Centrifuge *SET & LOCK STAT CENTRIFUGE*

Part # 76795-572

Discover a better way to spin STAT with the VWR® STAT 6 Clinical Centrifuge. Minimize lab error, reduce turnaround time, and streamline training with the VWR® STAT's simple controls and powerful performance.

Easy to Monitor

Prevent forgotten samples and reduce TAT with LED lid lighting. Bright lights indicate when the centrifuge is ready to load (off), running (on), or done (flashing). Lights continue flashing until the lid is opened.

Easy to Use

The simple two-button interface makes it easy to select a preset or program your custom validated STAT cycles. Use all three presets, or set and lock to a single cycle for maximum repeatability. Spin 75-100 mm (3-10 mL) tubes with the same rotor and universal tube holders.

Eliminate Wait Time

Further shorten TAT by using a centrifuge array in accessioning to eliminate queuing. The array takes advantage of the VWR® STAT 6's compact size to allow samples to be spun as soon as they arrive.

Safe and Durable

Enjoy years of worry-free operation with the VWR® STAT 6's brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.

A Great Value


Each centrifuge comes ready to run with a horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.




2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

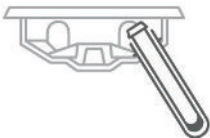

Rotation Angle

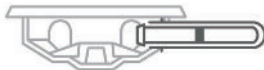

Applications	Max Speed	Max Capacity	Included	Warranty
STAT Chemistry, Coag, PPP	4,000 xg Max RCF 5,300 Max RPM	6 x 75-100 mm 3-10 mL	Horizontal Rotor, Tube Holders	2 years Rotor: lifetime warranty

Designed, built, and supported in the USA

VWR® STAT 6 Clinical Centrifuge



Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	76795-572
Replacement 75-100 mm Tube Holders/Buckets	03-1-0007-0046K (6 qty)
Adapters, 0.5 to 1 mL	7713068K (6 qty)
Adapters, 1.5 to 2 mL	7713065K (6 qty)
1 yr Warranty Extension	03-0-0007-0003
2 yr Warranty Extension	03-0-0007-0004
3 yr Warranty Extension	03-0-0007-0005

Factory Settings			
	RPM	Time	G-Force
1st Setting/Cycle	5,300	3	4,000
2nd Setting/Cycle	5,300	5	4,000
3rd Setting/Cycle	3,800	7	2,000

Specifications	
Sound level	61 dBA
Dimensions (W x D x H)	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	12 lbs (5.4 kg)
Power	95 to 253 VAC – 50/60 Hz
Max Speed/G-Force	5300 RPM/4000 xg



vwr.com | sales@vwr.com | 1-800-932-5000

VWR® STAT 12 Clinical Centrifuge *PROGRAMMABLE STAT CENTRIFUGE*

Part # 76795-576

Always know when your STAT samples are ready with the VWR® STAT 12 Clinical Centrifuge. The unique flashing lid lights (patent pending) make it easy to track cycle status and keep TAT down. Customize settings and check every detail on the digital display, and save up to 10 custom cycles to produce high-quality STAT samples.



Easy to Monitor

LED lid lighting indicates when the centrifuge is ready (off), running (on), or done (flashing), reducing forgotten samples and keeping your average TAT down. The VWR® STAT's digital display shows cycle parameters and time remaining.

Easy to Use

Set and save up to 10 programmable time, speed, and g-force settings using the user-friendly interface. The VWR® STAT 12 Clinical Centrifuge arrives pre-programmed with the most common STAT settings and can spin 75-100 mm (3-10 mL) tubes with the same rotor and universal tube holders.

Eliminate Wait Time

The VWR® STAT 12 Clinical Centrifuge is so compact that it can be used in an array in accessioning. The array will allow spinning samples as soon as they arrive, further shortening your TAT. For a more compact unit, consider the 6-tube VWR® STAT 6 Clinical Centrifuge

Safe and Durable

Enjoy years of worry-free operation with the STAT brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each VWR® STAT 12 centrifuge comes fully equipped and ready to run with a horizontal rotor, universal tube holders, and 2-year warranty included in the affordable price.


2-Stage Horizontal Rotor

Our patented swing-out horizontal rotor design incorporates a unique test tube holder that produces the ideal horizontally separated samples.

Loading Angle

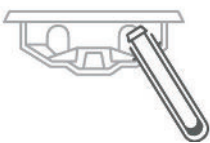

Rotation Angle

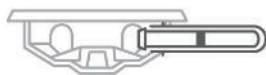

Applications	Max Speed	Max Capacity	Included	Warranty
STAT Chemistry, Coag, PPP	4,000 xg Max RCF 5,200 Max RPM	12 x 75-100 mm 3-10 mL	Horizontal Rotor, Tube Holders	2 years Rotor: lifetime warranty

VWR® STAT 12 Clinical Centrifuge

LED lid lights indicate cycle status: ready, running, done

QR code for easy access to support videos, manuals, and more



View RPM and remaining spin time on the digital display

Program and name your custom validated cycles

Store up to 10 settings

Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	76795-576
Replacement 75-100 mm Tube Holders/Buckets	03-1-0007-0046K (6 qty)
Adapters, 0.5 to 1 mL	7713068K (6 qty)
Adapters, 1.5 to 2 mL	7713065K (6 qty)
1 yr Warranty Extension	03-0-0007-0003
2 yr Warranty Extension	03-0-0007-0004
3 yr Warranty Extension	03-0-0007-0005

Factory Settings			
	RPM	Time	G-Force
1st Setting/Cycle	5,200	3 min	4,000 xg
2nd Setting/Cycle	5,200	5 min	4,000 xg
3rd Setting/Cycle	3,600	7 min	2,000 xg

Specifications	
Sound level	61 dBA
Dimensions (W x D x H)	13 x 15 x 9 in 33 x 38 x 23 cm
Weight	30 lbs (13.7 kg)
Power	95 to 253 VAC – 50/60 Hz
Max Speed/G-Force	5200 RPM/4000 xg



vwr.com | sales@vwr.com | 1-800-932-5000

VWR® STAT Coag Clinical Centrifuge

SET-AND-LOCK STAT COAG CENTRIFUGE

Part # 76795-588

Reduce turnaround time (TAT) and produce high-quality coags in minutes with the VWR® STAT Coag Clinical Centrifuge. The user-friendly set-and-lock controls streamline STAT lab operations and eliminate errors. Simply select your desired runtime from the STAT Coag's three preset cycles and spin with the press of a button.



Easy to Monitor

Always know when samples are ready with the VWR® STAT's LED lid lighting system (patent pending). Bright lights indicate when the centrifuge is ready to load (off), running (on), or done (flashing) and continue flashing until the lid is opened.

Easy to Use

Stop wasting time adjusting settings or navigating complicated menus. The VWR® STAT's simple interface makes it easy to select a preset or program your custom validated cycles with the press of a button. Spin tubes up to 75 mm (4 mL) with the same rotor and universal tube holders.

Eliminate Wait Time

Further shorten TAT by using a centrifuge array in accessioning. The array takes advantage of the VWR® STAT's compact size to eliminate queuing. Its light weight and small footprint allow it to be easily moved throughout the lab to address bottlenecks.

Safe and Durable

Enjoy years of worry-free operation with the brushless motor and rugged, reinforced carbon fiber components. Enhance lab safety with imbalance protection and a clear, shatterproof lid, which locks until the centrifuge has stopped completely.


A Great Value

Each centrifuge comes ready to run with fixed angle rotor and 2-year warranty included in the affordable price.


Fixed Angle Rotor

The rotor remains at a 45° angle during processing for optimal PPP formation.

Loading Angle

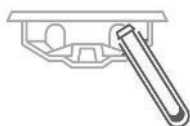

Rotation Angle

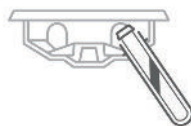

Applications	Max Speed	Max Capacity	Included	Warranty
STAT Coagulation, Platelet Poor Plasma (PPP)	4,400 xg Max RCF 6,600 Max RPM	12 x up to 75 mm 4 mL	Fixed Angle Rotor	2 years Rotor: lifetime warranty

Designed, built, and supported in the USA

VWR® STAT Coag Clinical Centrifuge



Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	76795-588
1 yr Warranty Extension	03-0-0007-0003
2 yr Warranty Extension	03-0-0007-0004
3 yr Warranty Extension	03-0-0007-0005

Factory Settings			
	RPM	Time	G-Force
1st Setting/Cycle	6,600	2	4,400
2nd Setting/Cycle	6,600	3	4,400
3rd Setting/Cycle	5,200	5	2,700

Specifications	
Sound level	61 dBA
Dimensions (W x D x H)	12 x 14 x 9 in 30 x 36 x 23 cm
Weight	12 lbs (5.4 kg)
Power	95 to 253 VAC – 50/60 Hz
Max Speed/G-Force	6600 RPM/4400 xg



vwr.com | sales@vwr.com | 1-800-932-5000

VWR SEROLOGY

VWR® SEROLOGY 12 Clinical Centrifuge **BLOOD BANKING CENTRIFUGE**

Part # 76795-586

The VWR® SEROLOGY 12 Clinical Centrifuge is carefully designed to facilitate high-quality blood grouping, typing, cross matching, and other cell-washing procedures for a wide range of tube sizes. Based on the world's most popular blood banking centrifuge and using the same rotors, the VWR® SEROLOGY 12 Clinical Centrifuge ensures samples are efficiently processed, with no individual handling necessary.

Easy to Use

The fully digital display makes it easy to customize cycle parameters in 5 or 15 second increments. Save up to 10 cycles, enabling easy spinning and simplifying validation. Spin 6 capped or 12 uncapped 75 mm blood tubes with minimal handling (up to 100 mm tubes with additional rotor). Tube positions are numbered for quick identification.

Easy to Monitor

See cycle status at a glance with the VWR® SEROLOGY 12 Clinical Centrifuge's LED lid lighting (patent pending). Bright lights indicate when the centrifuge is ready to load (off), running (on), or done (flashing). Lights continue flashing until the lid is opened.

No Individual Handling

Tubes can be batch filled, centrifuged, transported, incubated, and emptied using the interchangeable heads that also function as holders.

Safe and Durable

Enhance lab safety with imbalance protection and a shatterproof lid, which locks until the rotor has stopped completely. The brushless motor and reinforced carbon fiber components ensure years of operation with no routine maintenance.

A Great Value

Unbeatable quality for any budget. The VWR® SEROLOGY 12 Clinical Centrifuge is ready to run with a 2-stage swing-out rotor and 2-year warranty included in the affordable price.



Uses Sero-Fuge™ Rotors

The VWR® SEROLOGY 12 Clinical Centrifuge is compatible with all rotors from the world's most popular blood banking centrifuge.



Sold separately



Included with purchase

Applications	Max Speed	Max Capacity	Included	Warranty
Cell Washing Blood Grouping & Typing Cross Matching	1,040 xg Max RCF 3,600 Max RPM	6 capped or 12 uncapped 75mm tubes*	12 tube rotor <i>*6 place metal rotor (sold separately) required for 100 mm tubes</i>	2 years

VWR® SEROLOGY 12 CLINICAL CENTRIFUGE



QR code for easy access to support videos, manuals, and more

Easily verify spin settings with the clear digital display

Adjust time in 5 second blocks

Airflow cooling maintains temperature during processing



Rotor Compatibility

Transition to a new centrifuge without interrupting workflow. The VWR® SEROLOGY 12 Clinical Centrifuge is compatible with all rotors from the world's most popular blood banking centrifuge.

Free 30-Day Evaluation (U.S. Only)

Part Numbers	
Centrifuge	76795-586

Specifications	
Sound level	61 dBA
Dimensions (W x D x H)	13 x 15 x 9 in (33 x 38 x 23 cm)
Weight	27 lbs (12.3 kg)
Maximum Speed	3600 RPM
Maximum RCF	1040 xg
Voltage	95-253 VAC
Frequency	50/60 Hz



6 place metal rotor
0526-001-000 (1 qty)
for 100mm tubes
(sold separately)



12 place rotor
02-001-0-0019 (1 qty)
for 75mm tubes
(1 included)



vwr.com | sales@vwr.com | 1-800-932-5000