

# HORIZON

## DRUCKER HORIZON 6 FA *FIXED ANGLE VETERINARY CENTRIFUGE*



The HORIZON 6 FA is a compact fixed angle veterinary centrifuge. Modeled after the world's best-selling 642E clinical laboratory centrifuge, the HORIZON 6 FA is easy to use, with simple controls designed for maximum repeatability and a single preset for all common veterinary needs. Its small footprint makes it an easy fit in even the tightest labs.

### Easy to Use

Spin samples with the press of just one button. The HORIZON 6 FA processes 125 mm tubes and features three presets for all common veterinary needs.

### Built to Last

No routine maintenance. The brushless motor and unique convection-based cooling system ensure years of high-quality separations.

### Engineered for Safety

Protect your lab with the HORIZON 6 FA clear, safety locking lid. The unique convection-based cooling system and brushless motor ensure years of worry-free operation.


### A Great Value

Get the HORIZON 6 FA at a price that fits any budget. The unit arrives ready for use with rotor, tube holders, and 2-year warranty included in the affordable price.


### Fixed-Angle Rotor Included

The rotor remains at a 45° angle during processing. Spin 75 mm, 100 mm, and 125 mm tubes in specially designed tube holders.

#### Loading Angle

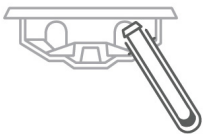

#### Rotation Angle

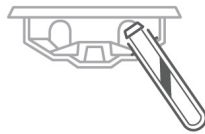

Tubes	Applications	Centrifugation	Included	Warranty
6 15 mL/125 mm	Veterinary	Fixed Angle 1850 xg Max RCF 3900 Max RPM	Rotor Tube Holders	2 years Rotor: lifetime warranty Extended warranty available

Designed, built, and supported in the USA

Scan to buy:



# HORIZON 6 FA

Safety lid automatically locks during operation and will not open until the centrifuge has stopped

Lights illuminate to indicate running, locked, unlocked

Simple two-button controls



Reinforced guard bowl for maximum safety

Brushless motor ensures quiet, maintenance-free operation

## Free 30-Day Evaluation (U.S. Only)

Standard Preset Cycles			
Cycle 1	2 min	2000 RPM	500 xg
Cycle 2	10 min	3900 RPM	1600 xg
Cycle 3	6 min	3900 RPM	1550 xg

Part Numbers	
Centrifuge	00-276-223-000
125 mm tube holder (1 qty)	7713032
0.5-1 mL Purple Adapters (6 qty)	03-1-0007-0059K
Microhematocrit Insert Kit	02-004-0-0029K
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Additional Accessories for Purchase

Specifications	
Sound level	64 dBA
Dimensions	12 x 14 x 9 in 30 x 35.5 x 23 cm
Weight	12 lbs (5.4 kg)
Max Acceleration	10 seconds
Max Deceleration	20 seconds
Power	100 to 240 VAC – 50/60 Hz



druckerdiagnostics.com  
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)  
+1-814-692-7661