

HORIZON

DRUCKER HORIZON 6 FA *FIXED ANGLE VETERINARY CENTRIFUGE*

Part # 00-276-226-000



The HORIZON 6 FA is a compact fixed angle veterinary centrifuge. Modeled after the world's best-selling 642E clinical laboratory centrifuge, the HORIZON 6 FA is easy to use, with simple controls designed for maximum repeatability and a single preset for all common veterinary needs. Its small footprint makes it an easy fit in even the tightest labs.

Easy to Use

Spin samples with the press of just one button. The HORIZON 6 FA processes tubes up to 16 x 125 mm in the included orange carriers and features three presets for all common veterinary needs.

Built to Last

No routine maintenance. Brushless motor and unique convection-based cooling system ensures years of high-quality separations.

Engineered for Safety

Protect your lab with imbalance detection and a shatterproof lid which can only be opened when the rotor is stopped, even if the unit loses power.


A Great Value

Unbeatable quality for any budget. Each HORIZON 6 FA centrifuge is ready to run with fixed angle rotor, black carriers, and 2-year warranty included in the affordable price.


Fixed Angle Rotor Included

The rotor remains at a 45° angle during processing. Spin 16 x 125 mm tubes in the specifically designed carriers.

Loading Angle

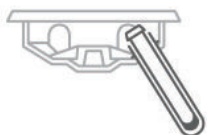

Rotation Angle

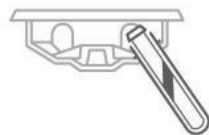

Applications	Max Speed	Max Capacity	Included	Warranty
Veterinary	3900 Max RPM 1850 xg Max RCF	6 x up to 16 x 125 mm or 15 mL tubes	Fixed Angle Rotor 6 x 16 x 125 mm Black Carrier	2 years Extended warranty options available

Designed, built, and supported in the USA

HORIZON 6 FA

Scan to buy:



Safety lid automatically locks during operation and will not open until the centrifuge has stopped

Lights illuminate to indicate running, locked, unlocked

Simple two-button controls



Reinforced guard bowl for maximum safety

Brushless motor ensures quiet, maintenance free operation



not included with centrifuge purchase

Spin Microhematocrit Tubes and 1-1.5 mL Microtainer Tubes

Spin up to 8 microhematocrit tubes with out microhematocrit adapter kit, designed to fit either 40mm or 75mm microhematocrit tubes. The adapter kit includes 2 orange microhematocrit inserts, 6 blue carries, 6 microtainer adapters and veterinary microhematocrit card reader.

Standard Preset Cycles	Time	RPM	RCF
Cycle 1	2 min	2000 RPM	500 xg
Cycle 2	10 min	3900 RPM	1600 xg
Cycle 3	6 min	3900 RPM	1550 xg

Part Numbers	
Centrifuge	00-276-226-000
16 x 125 mm Carrier, Black (6 qty)	7713032K
Microhematocrit Adapter Kit	02-004-0-0029K
1 yr Warranty Extension	03-0-0007-0001
2 yr Warranty Extension	03-0-0007-0002
3 yr Warranty Extension	03-0-0007-0006

Specifications	
Sound Level	64 dBA
Max RPM	3900 RPM
Max G-Force	1850 xg
Dimensions (W x D x H)	12 x 14 x 9 in 30.5 x 35.6 x 22.9 cm
Weight	12 lbs (5.4 kg)
Power	100 - 240 (+/- 10%) – 50/60 Hz



druckerdiagnostics.com
sales@druckerdiagnostics.com

+1-866-265-1486 (U.S. only)
+1-814-692-7661